


Analysis of Benzene Derivatives

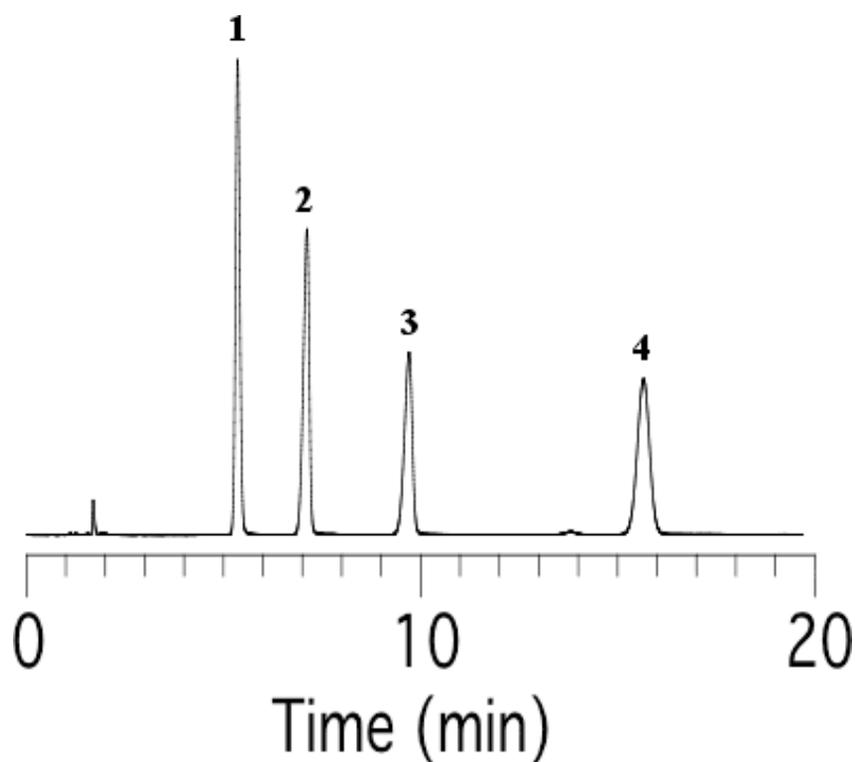

Column : Inertsil Ph-3
 5 μ m
 150 x 4.6 mm I.D.

Column Cat. No. : 5020-01302

Eluent : A) CH₃CN
 B) H₂O
 A/B = 30/70, w/w

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :

Solutes : 1. Acetophenone
 2. Benzene
 3. Toluene
 4. Naphthalene