


# Analysis of Aromatic Compounds

Data No. LA437-0000

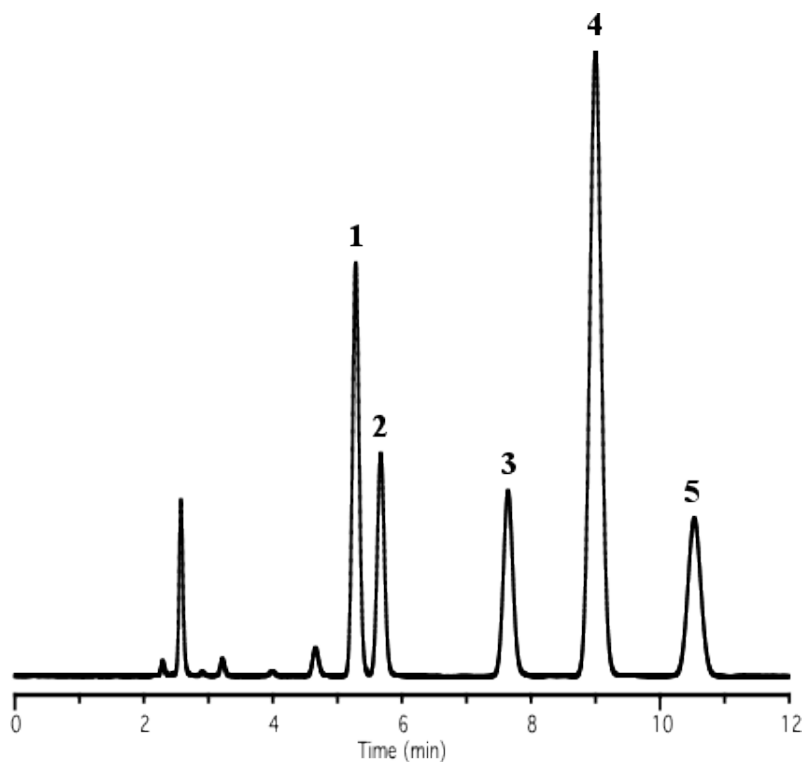

Column	: Inertsil ODS-3 5 µm 150 x 4.6 mm I.D.	Solutes	: 1. Naphthalene 2. Ethylbenzene 3. n-Propylbenzene 4. Anthracene 5. n-Butylbenzene
Column Cat. No.	: 5020-01731		
Eluent	: A) CH <sub>3</sub> CN B) H <sub>2</sub> O A/B = 70/30, v/v		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 254 nm		
Injection Vol.	:		
Sample	:		