


## Analysis of Pyridine

Data No. LA450-0000

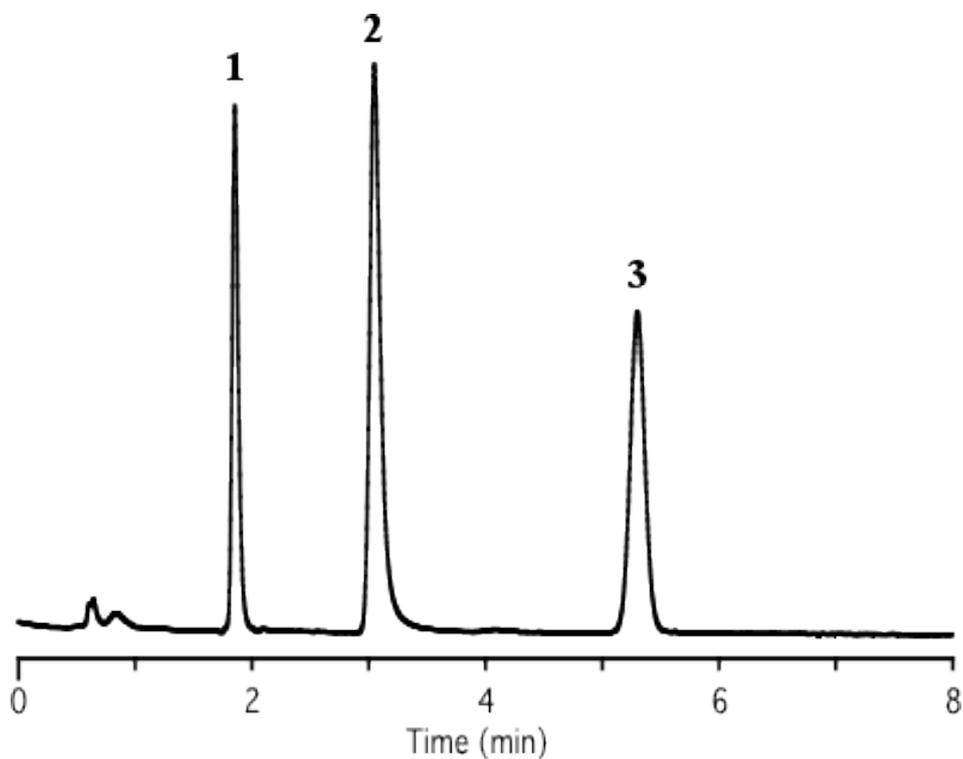

Column	: Inertsil CN-3 5 $\mu$ m 150 x 4.6 mm I.D.	Solutes	: 1. Uracil 2. Pyridine 3. Phenol
Column Cat. No.	: 5020-01940		
Eluent	: A) CH <sub>3</sub> OH B) H <sub>2</sub> O A/B = 30/70, v/v		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 254 nm		
Injection Vol.	:		
Sample	:		