


## Analysis of Benzene Derivatives

Data No. LA452-0000

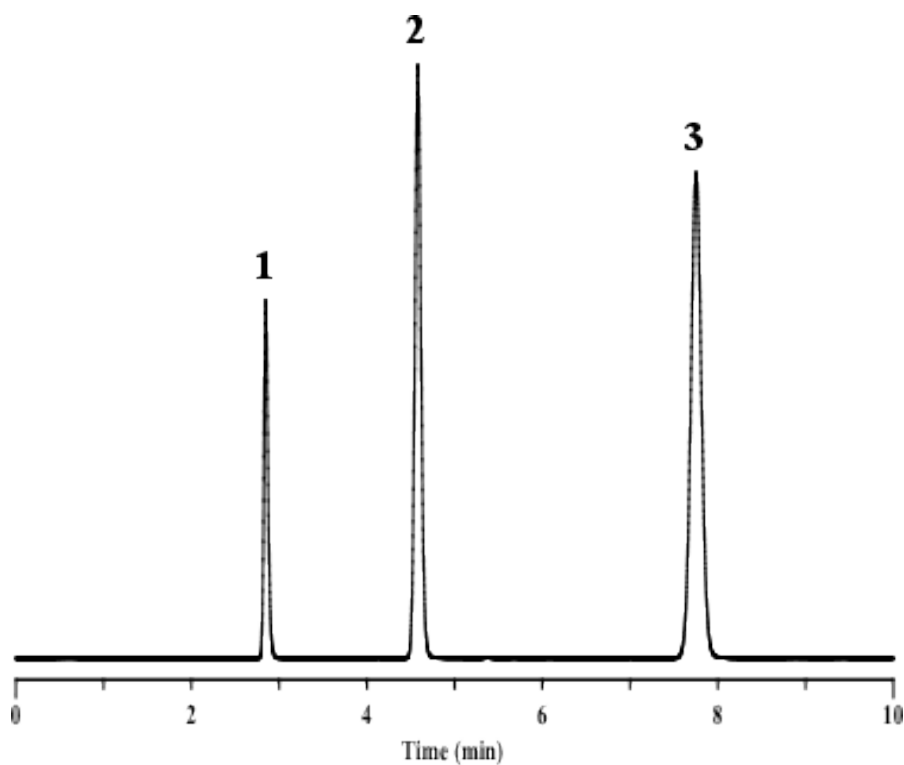

Column	: Inertsil Diol 5 µm 250 x 4.6 mm I.D.	Solutes	: 1. Benzene 2. Nitrobenzene 3. o-Nitroanisoie
Column Cat. No.	: 5020-05646		
Eluent	: A) n-Hexane B) Ethanol A/B = 98/2, v/v		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 254 nm		
Injection Vol.	:		
Sample	:		