


Analysis of Aromatic Compounds

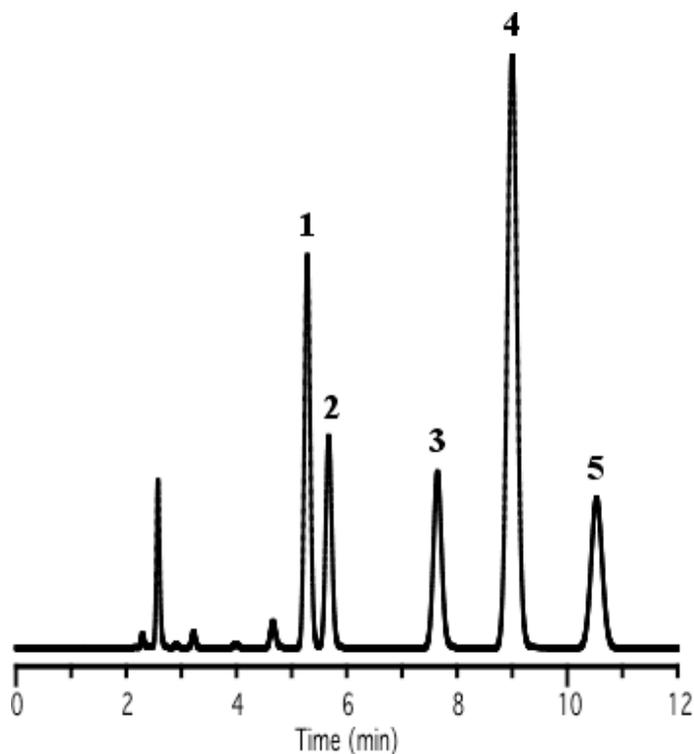

Column : Inertsil ODS-3
 5 μ m
 150 x 4.6 mm I.D.

Column Cat. No. : 5020-01731

Eluent : A) CH₃CN
 B) H₂O
 A/B = 70/30, v/v

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :

Solutes : 1. Naphthalene
 2. Ethylbenzene
 3. n-Propylbenzene
 4. Anthracene
 5. n-Butylbenzene