


Analysis of Aromatic Compounds

Data No. LA462-0000

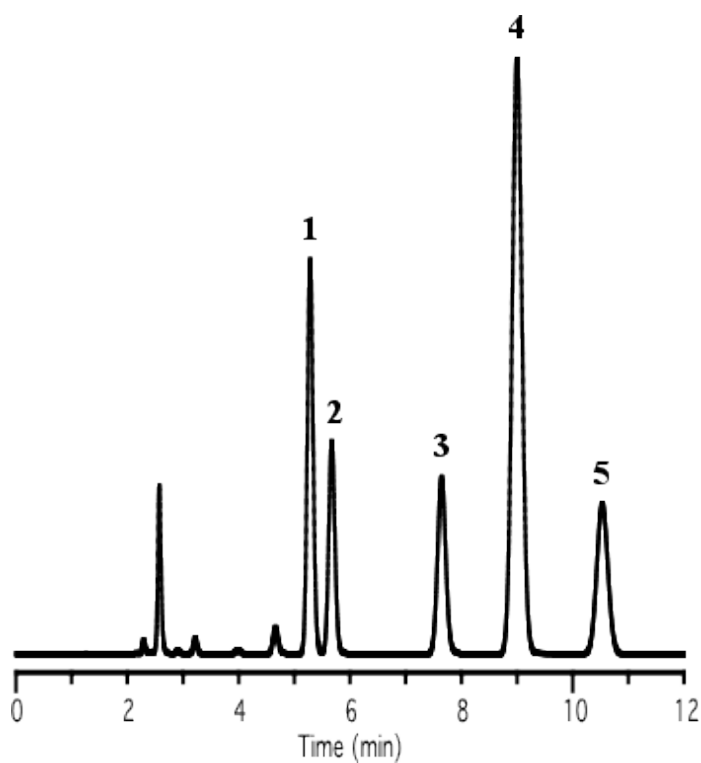

Column : Inertsil ODS-3
5 µm
150 x 4.6 mm I.D.

Column Cat. No. : 5020-01731

Eluent : A) CH₃CN
B) H₂O
A/B = 70/30, v/v

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :

Solutes : 1. Naphthalene
2. Ethylbenzene
3. n-Propylbenzene
4. Anthracene
5. n-Butylbenzene