


InertSearch™ for LC

Inertsil® Applications

Analysis of 5-Fluorouracil

Data No. LA548-0696

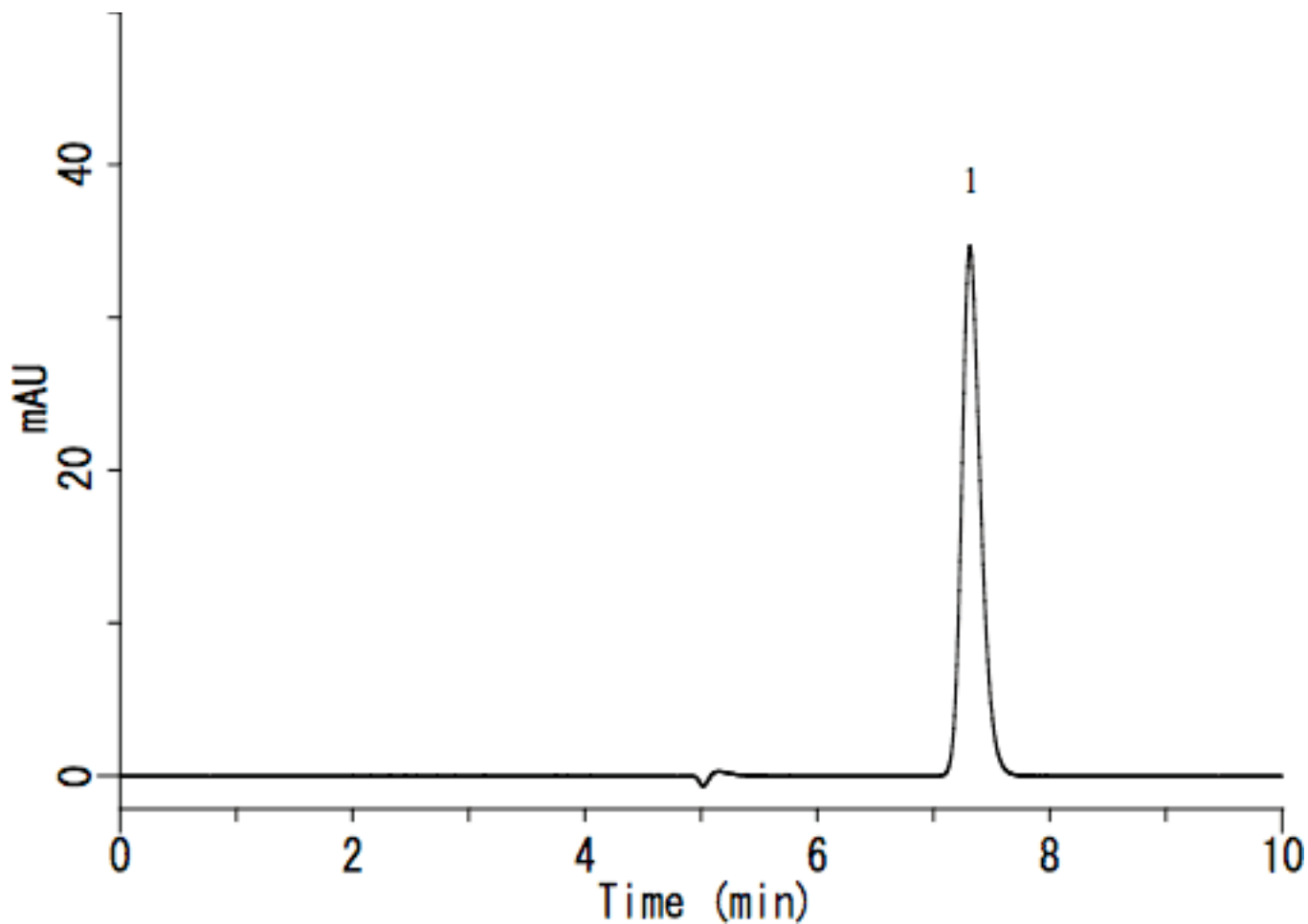

System : GL-7400 HPLC system Solutes : 1. 5-Fluorouracil(10 mg/L)
Column : Inertsil ODS-3
 5 μ m,
 250 \times 3.0 mm I.D.
Column Cat. No. :
Eluent : A) CH₃CN
 B) H₂O
 A/B = 2/98, v/v (Mixed by a gradient mixer)
Flow rate : 0.4 mL/min.
Col. Temp. : 40 $^{\circ}$ C
Detection : PDA 265 nm
Injection Vol. : 5 μ L
Sample : Anti cancer drug