
Analysis of Biodiesel

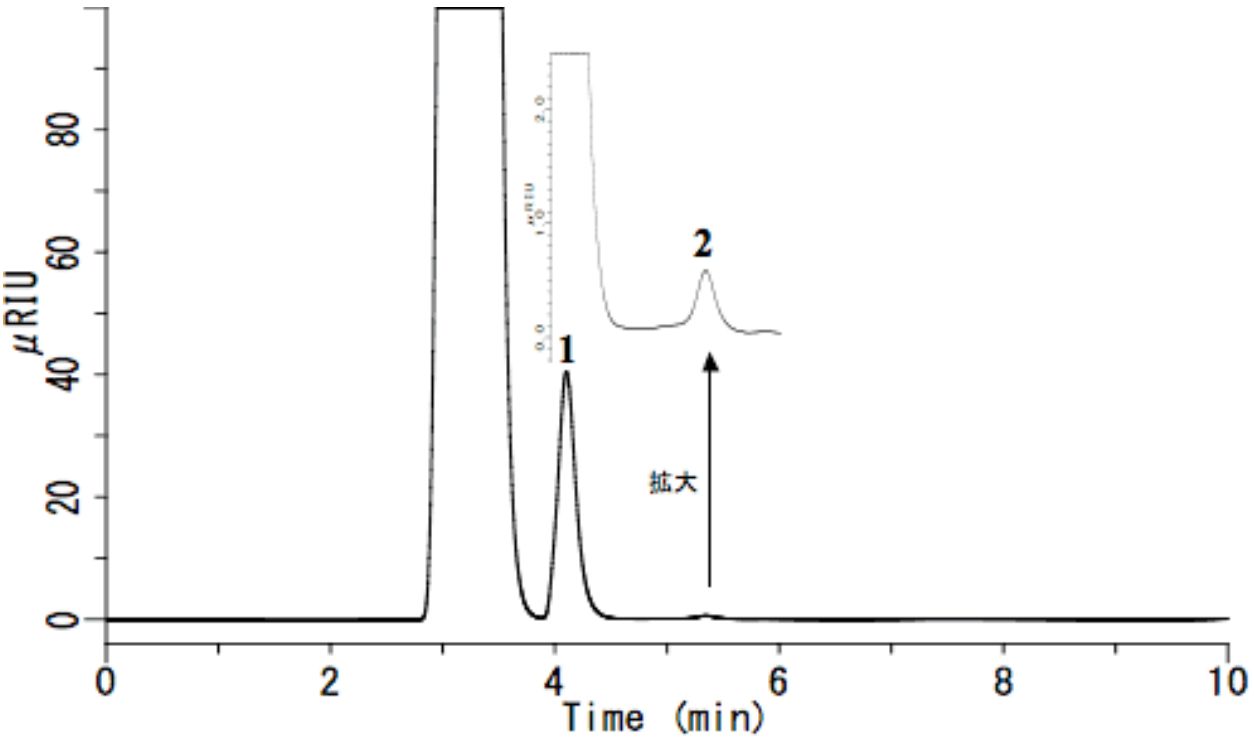

System	:GL-7400 system	Solutes	:1. Methyl Stearate (2000 mg/L)
Column	:Inertsil SIL-100A (5µm, 250×4.6 mm I.D.)		2. Trilinolein(20 mg/L)
Column Cat. No.:			
Eluent	:A) n-Hexane		
	B) 2-Propanol		
	A/B = 996/4, v/v		
Flow rate	:1.0 mL/min		
Col. Temp.	:40 °C		
Detection	:RI (GL-7454 RI detector)		
Injection Vol.	:10 µL		
Sample	:Methyl Stearate		
	Trilinolein		