


## Analysis of Tryptic Digest of Cytochrom C

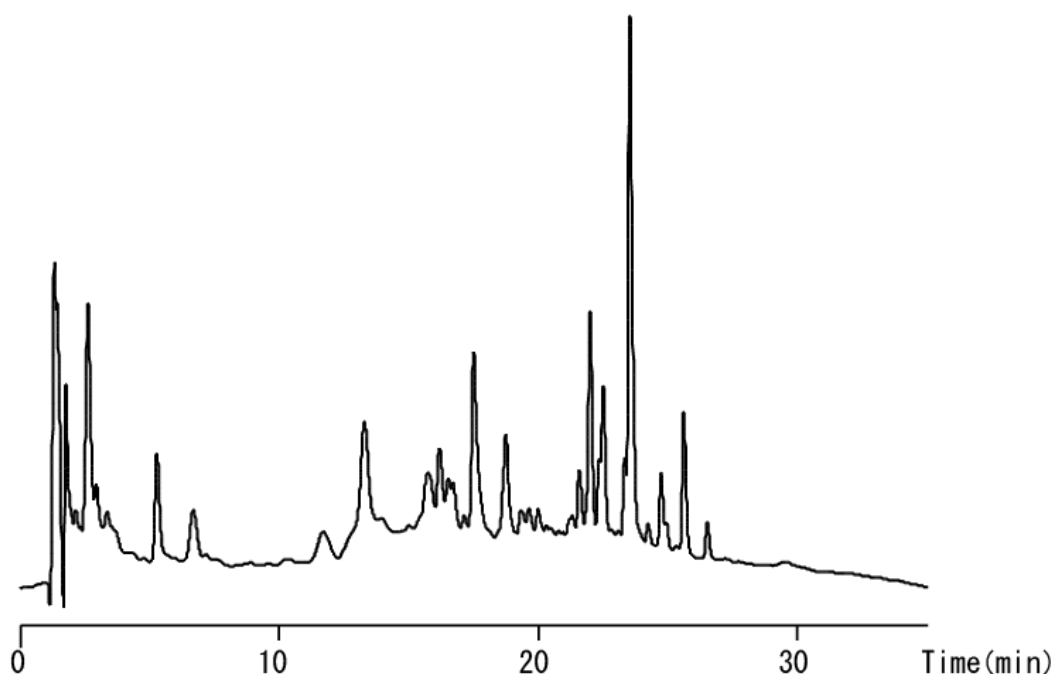

### Conditions

**Column** : Inertsil Peptides C18 (4  $\mu$ m, 150 x 2.1 mm I.D.)  
**Column Cat. No.** : 5020-08025  
**Eluent** : A) 0.1 % TFA in (CH<sub>3</sub>CN/H<sub>2</sub>O = 90/10, v/v)  
               B) 0.1 % H<sub>2</sub>O  
               A/B = 10/90 – 10 min – 10/90 – 40 min – 60/40, v/v  
**Flow Rate** : 0.3 mL/min  
**Col. Temp.** : 30 °C  
**Detection** : UV 210 nm  
**Injection Vol.** : 2  $\mu$ L  
**Sample** : Tryptic Digest of Cytochrom C