


InertSearch for LC Data No. LA672-0699

Analysis of Protein (Inertsil WP300 Diol)

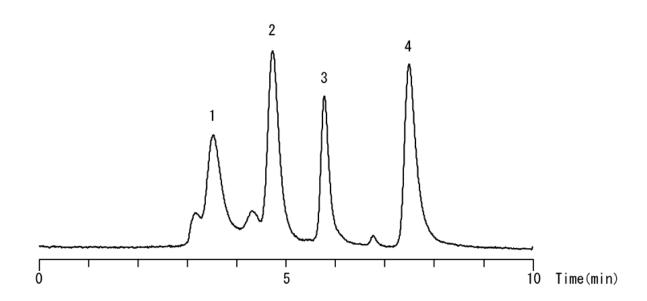

Conditions

Column : Inertsil WP300 Diol (5 μm, 250 x 4.6 mm I.D.)

Column Cat. No.: 5020-05746

Eluent : $0.1 \text{ M Na}_2\text{HPO}_4 \text{ (pH 6.9, NaH}_2\text{PO}_4\text{)}$

Analyte

1. Thyrogloblin	(0.25 mg/mL)
2. BSA	(0.25 mg/mL)
3. Insulin Chain A	(0.25 mg/mL)
4. Oxytocin	(0.21 mg/mL)