


InertSearch for LC Data No. LA730-0871

Analysis of Tranexamic acid

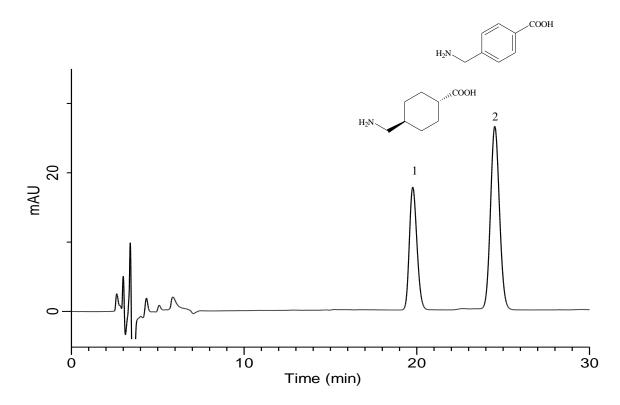

Conditions

System : GL-7400 HPLC system

Column : Inertsil ODS-3 (5μm, 250 x 4.6 mm I.D.)

Column Cat. No. : 5020-01732

Eluent : A) Phosphate buffer*

B) CH₃OH

A/B = 60/40, v/v

Flow Rate : 0.7 mL/min

Col. Temp. : 25 °C

Detection : UV 220 nm (GL-7450 UV detector)

Injection Vol. : 20 μL

* Phosphate buffer:

11.0 g of NaH_2PO_4 , 5 mL of triethyl amine and 1.4 g of sodium lauryl sulfate were dissolved in 500 mL of H_2O , and H_3PO_4 was added to the solution to adjust the pH value to 2.5. H_2O was added again to make up the solution to 600 mL.

Analyte

1. Tranexamic acid (1 mg/mL)
2. 4-(aminomethyl)benzoic acid (0.01 mg/mL)