


Analysis of Sugars (MonoCap Amide)

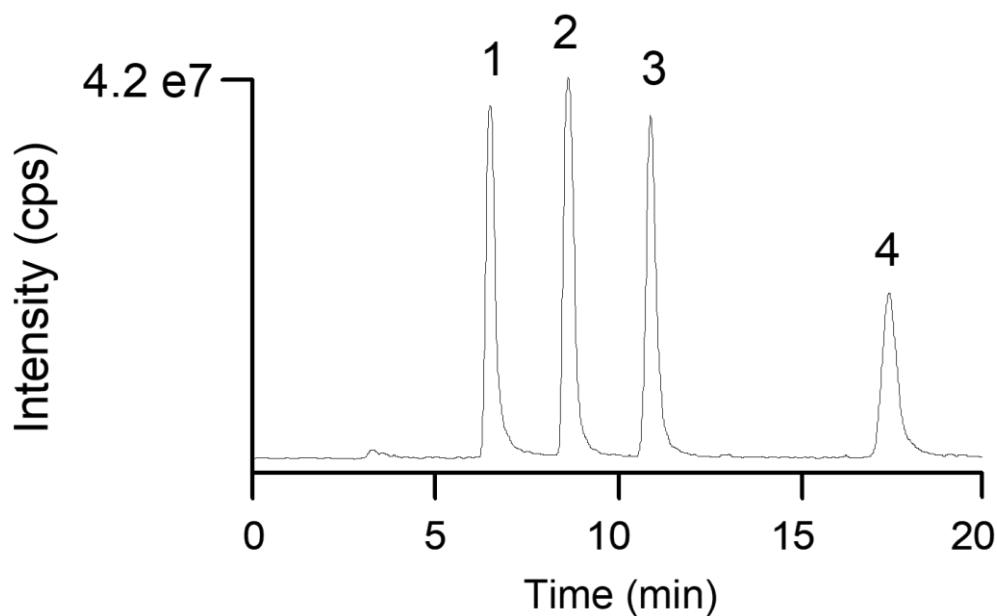

Conditions

System : AccuStream 700 series
Column : MonoCap Amide (250 x 0.2 mm I.D.)
Column Cat. No. : 5020-10172
Eluent : A) CH₃CN
 : B) 10 mM CH₃COONH₄
 : A/B = 80/20, v/v
Flow Rate : 2 µL/min
Col. Temp. : Ambient
Detection : Nano-ESI (Negative, 2000 V, MonoSpray FS 100 µm, TIC 180 – 530)
Injection Vol. : 0.01 µL

Analyte

1. Sorbitol	(20 µg/mL)
2. Sucrose	(20 µg/mL)
3. Trehalose	(20 µg/mL)
4. Raffinose	(20 µg/mL)