


InertSearch™ for LC

Inertsil® Applications

Analysis of Sugars (MonoCap Amide)

Data No. LA794-0721

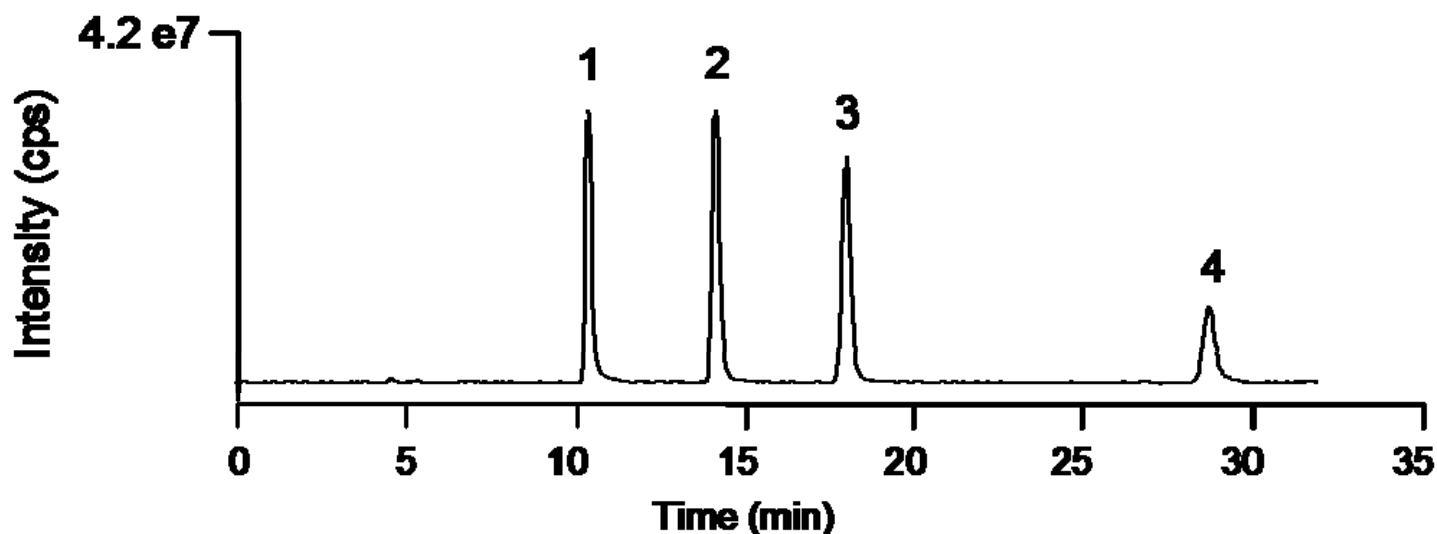

Conditions

System	: AccuStream 700 series	Analyte	
Column	: MonoCap Amide (500 x 0.2 mm I.D.)	1. Sorbitol	(20 µg/mL)
Column Cat. No.	: 5020-10173	2. Sucrose	(20 µg/mL)
Eluent	: A) CH ₃ CN B) 10 mM CH ₃ COONH ₄ A/B = 80/20, v/v	3. Trehalose	(20 µg/mL)
Flow Rate	: 2 µL/min	4. Raffinose	(20 µg/mL)
Col. Temp.	: Ambient		
Detection	: Nano-ESI (Negative, 2000 V, MonoSpray FS 100 µm, TIC 180 – 530)		
Injection Vol.	: 0.01 µL		