


## Analysis of Vitamin K (Inertsil ODS-4)

Data No. LA875-0894

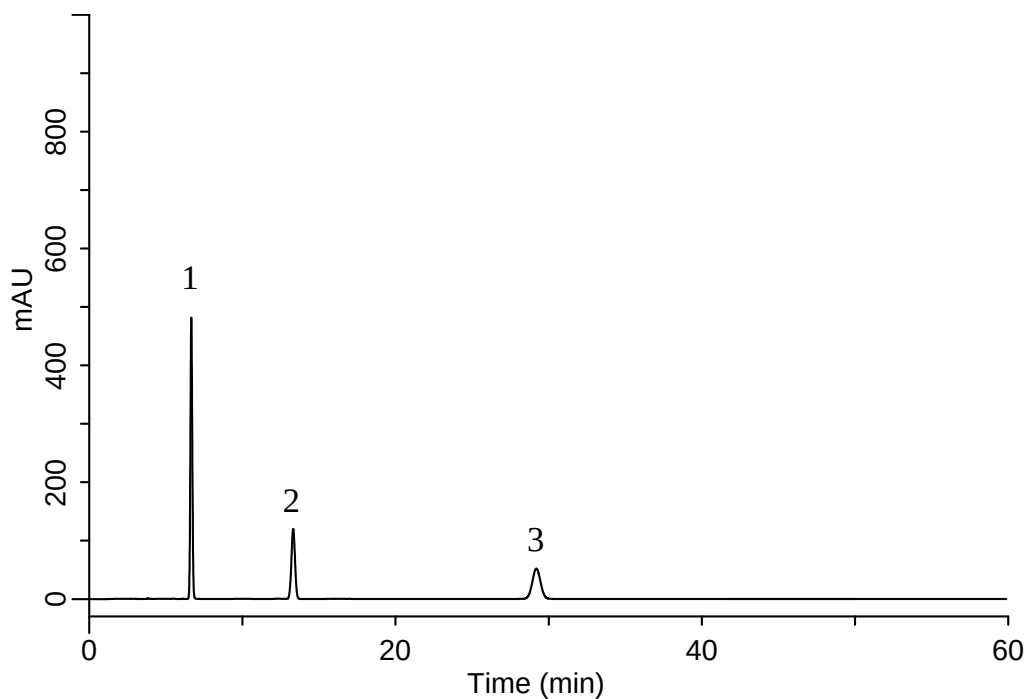

### Conditions

**System** : GL-7400 HPLC system  
**Column** : Inertsil ODS-4 (5  $\mu$  m, 150 x 4.6 mm I.D.)  
**Column Cat. No.** : 5020-03945  
**Eluent** : CH<sub>3</sub>CN  
**Flow rate** : 1.0 mL/min  
**Col. Temp.** : 40 °C  
**Detection** : UV 270 nm (GL-7452A PDA Detector)  
**Injection Vol.** : 5  $\mu$  L  
**Sample** : Vitamin K

### Solutes:

|    |                   |          |
|----|-------------------|----------|
| 1. | Vitamin K2 (MK-4) | 400 mg/L |
| 2. | Vitamin K1        | 200 mg/L |
| 3. | Vitamin K2 (MK-7) | 250 mg/L |