


Analysis of Sugar alcohols

Data No. LB008-0721

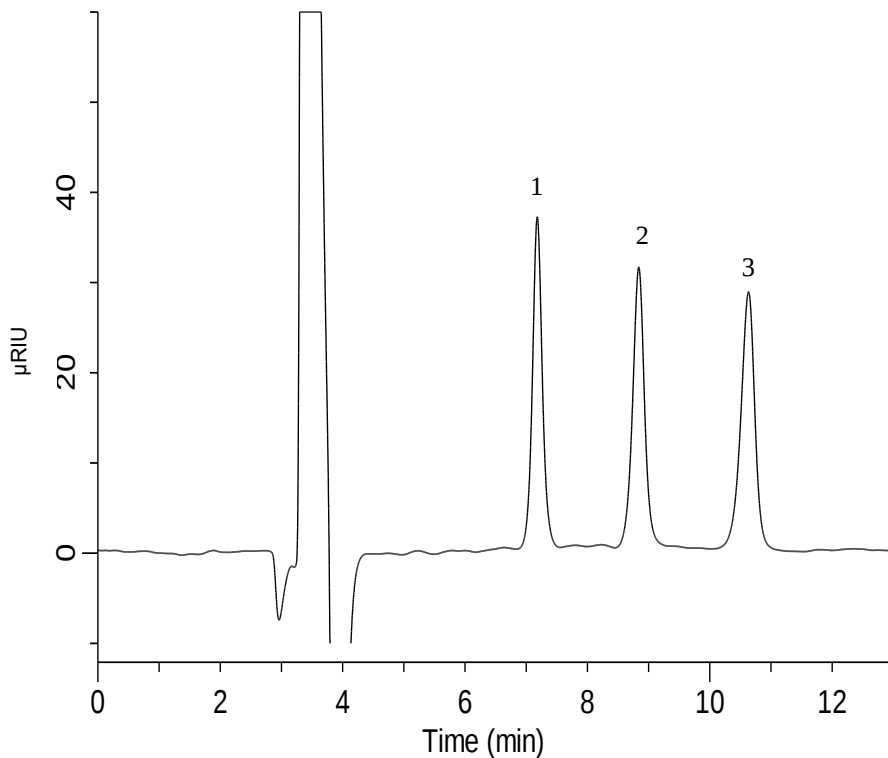

Conditions

System : GL-7400 HPLC system
Column : Inertsil Amide (5 μm, 250 x 4.6 mm I.D.)
Column Cat. No. : 5020-07836
Eluent : A) CH₃CN
 B) H₂O
 A/B = 75/25, v/v
Flow rate : 1 mL/min
Col. Temp. : 40 °C
Detection : RI (GL-7454 RI Detector)
Injection Vol. : 1 μL
Sample : Standard

Analyte:

1. meso-Erythritol	1 mg/mL
2. Xylitol	1 mg/mL
3. Sorbitol	1 mg/mL

