
Analysis of Linear alkylbenzene sulfonate

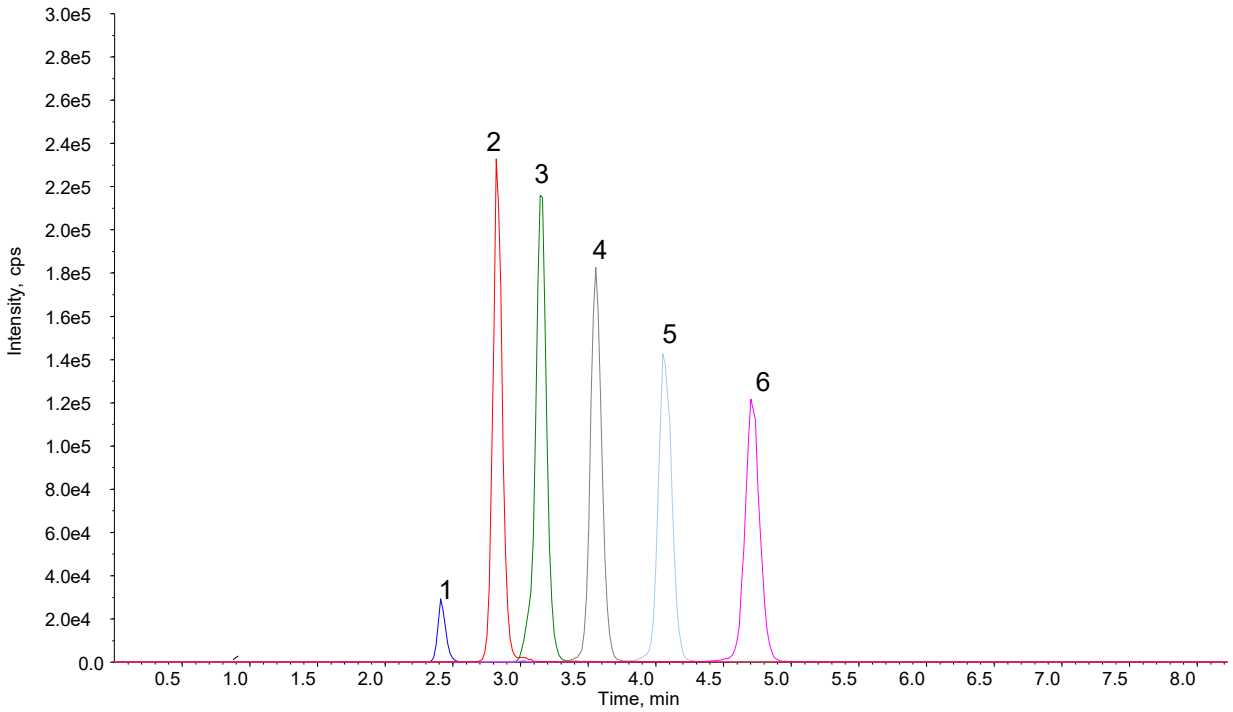

Conditions

System : LC800 HPLC System
4000 QTRAP®

Column : Inertsil C8-4 HP
(3 μm, 150×2.1 mm I.D.)

Column Cat. No. : 5020-03975

Eluent : A) CH_3CN
B) 0.1 % HCOOH + 50 mM $HCOONH_4$ in H_2O
A/B = 65/35, v/v

Flow Rate : 0.2 mL/min.

Col. Temp. : 40 °C

Detection : LC/MS/MS (4000 QTRAP® ESI, Negative, MRM)

CUR	CAD	IS	TEM	GS1	GS2
10	4	-4500	600	70	40

Injection Vol. : 5 μL

Analyte:

1. Sodium Octylbenzenesulfonate(C8) (I.S.)
2. Sodium Decylbenzenesulfonate(C10)
3. Sodium Undecylbenzenesulfonate(C11)
4. Sodium Dodecylbenzenesulfonate(C12)
5. Sodium Tridecylbenzenesulfonate(C13)
6. Sodium Tetradecylbenzenesulfonate(C14)

in CH_3CN/H_2O = 65/35, v/v
each 25 μg/L (I.S. 5 μg/L)