
Analysis of 2 kinds of Aldehydes (Detected by Pre-Column Method with DNPH)

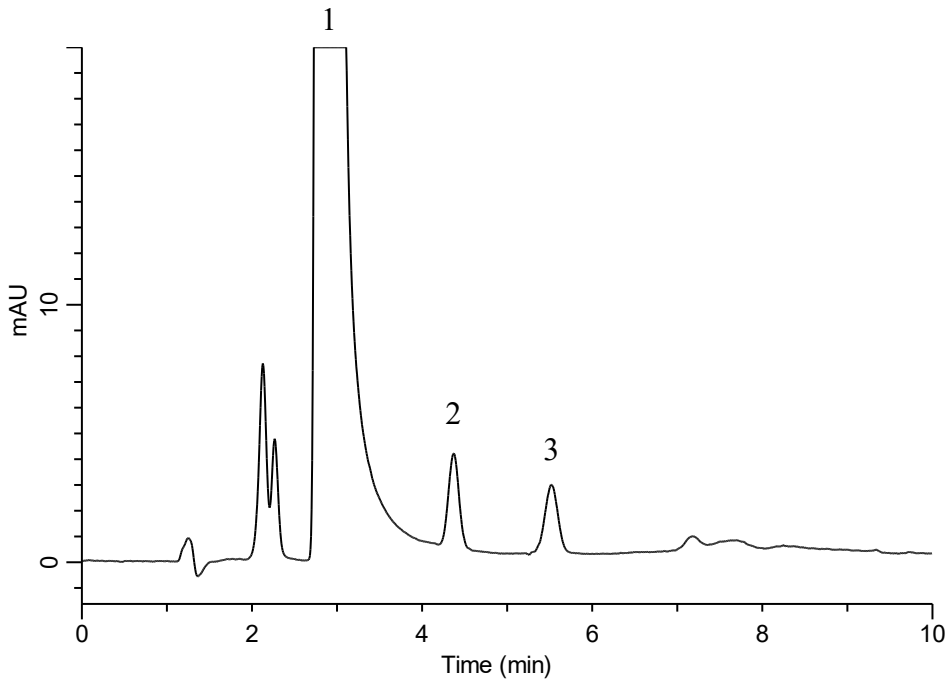

Conditions

System : GL-7400 HPLC system
Column : Inertsil ODS (5 μm, 150 x 4.6 mm I.D.)
Column Cat. No. : 5020-02124
Eluent : A) CH_3CN
B) H_2O
A/B = 60/40, v/v
Flow Rate : 1.0 mL/min
Col. Temp. : 40 ℃
Detection : UV 360 nm (GL-7452 PDA Detector)
Injection Vol. : 20 μL
Sample : Standard

Analyte:

1. DNPH
2. DNPH-Formaldehyde 50 μg/L
3. DNPH-Acetaldehyde 50 μg/L