


InertSearch for LC

Inertsil Applications

Analysis of Insulin

Data No. LB664-0812

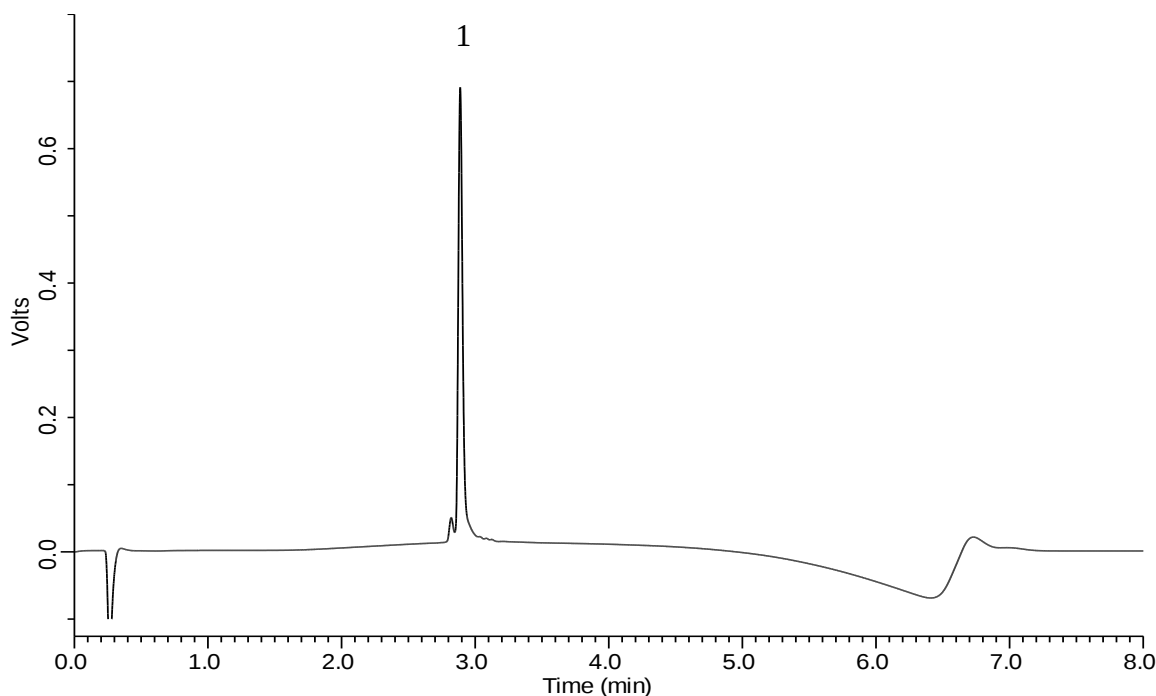

Conditions

Column Cat. No. : 5020-10818

Column : MonoSelect RP-mAb
(20 x 2.1 mm I.D.)

Eluent : A) 0.1%TFA in CH₃CN
B) 0.1%TFA in H₂O
A/B = 5/95-5 min-90/10-0.1 min-5/95-3 min

Flow Rate : 0.3 mL/min

Col. Temp. : 80 °C

Detection : UV 210 nm

Injection Vol. : 5 µ L

Sample : Standard

Analyte:

1. Insulin

0.1 mg/mL