


Analysis of D(-) -Mannitol

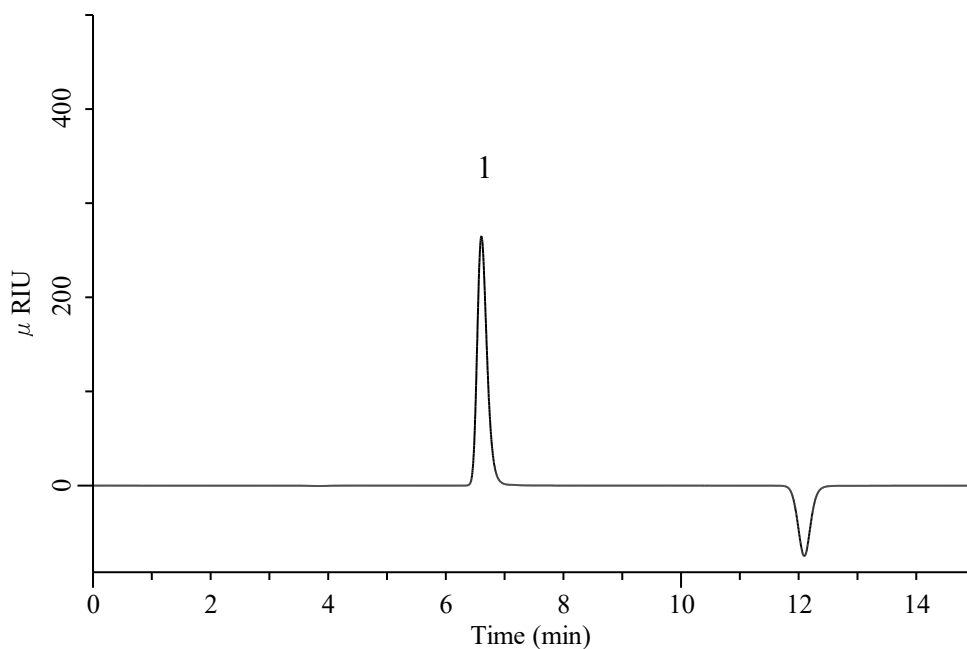

Conditions

System : Chromaster PLUS (HITACHI)
Column : InertSphere FA-1 (GL Sciences Inc.)
 (9 μm, 300 x 7.8 mm I.D.)
Column Cat. No. : 5020-11003
Eluent : H₂O
Flow Rate : 1.0 mL/min
Col. Temp. : 80 °C
Detection : RI
Injection Vol. : 80 μL
Sample : Standard in H₂O

Analyte:

1. D(-)-Mannitol 2.5 mg/mL

Tailing factor : 1.30 (≦ 2.0)

RSD of the peak area (%)(n=6) : 0.04 (≦ 2.0)