


Analysis of Absciscic acid

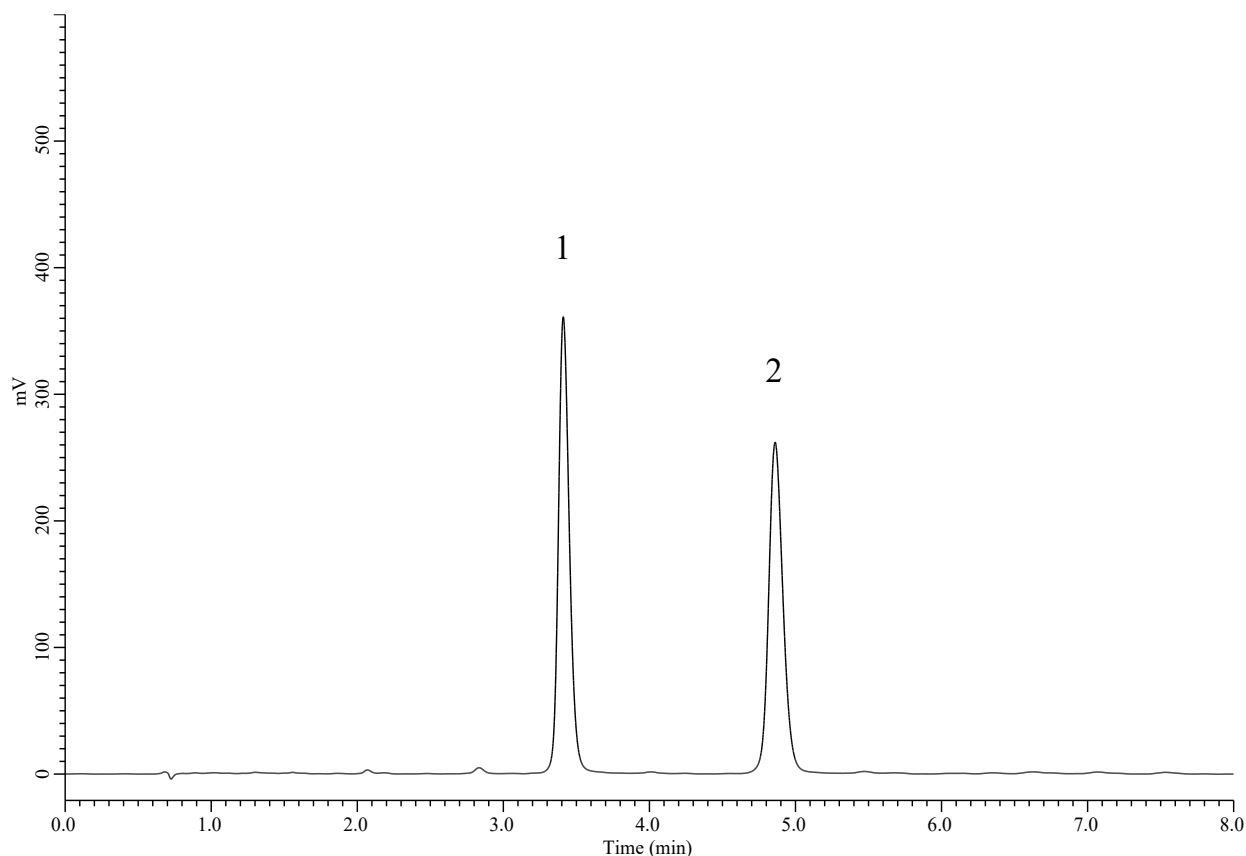

Conditions

Column : InertCore Biphenyl (GL Sciences Inc.)
(2.6 μm , 100 x 2.1 mm I.D.)

Column Cat. No. : 5020-17544

Eluent A) CH_3OH
B) 0.1% HCOOH in H_2O
: A/B = 40/60, v/v

Flow Rate : 0.3 mL/min

Col. Temp. : 40 $^{\circ}\text{C}$

Detection : UV 240 nm

Injection Vol. : 1 μL

Analyte:

1. *t-t*-Absciscic acid

2. *c-t*-Absciscic acid

(100 mg/L each)