


Analysis of Antiviral drugs

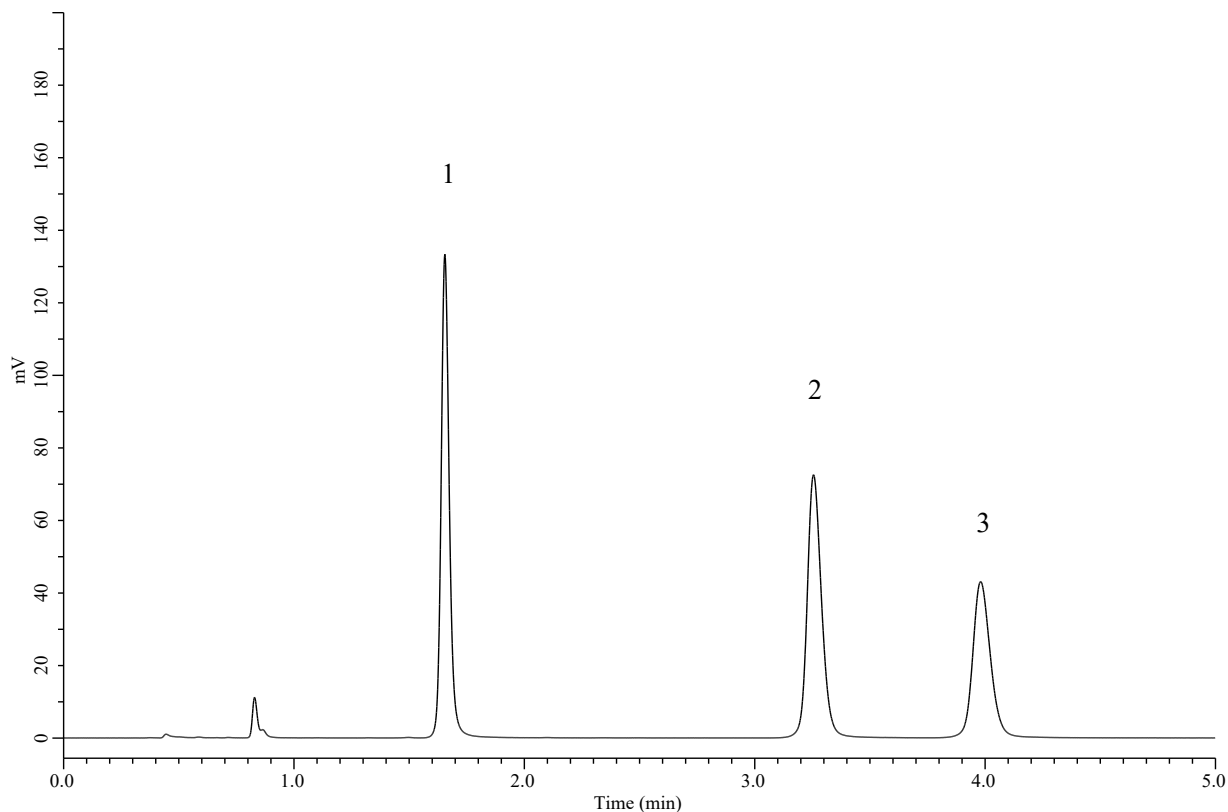

Conditions

Column : InertCore Biphenyl (GL Sciences Inc.)
 (2.6 μm , 100 x 2.1 mm I.D.)
Column Cat. No. : 5020-17544
Eluent : 10 mM HCOONH_4 in H_2O
Flow Rate : 0.4 mL/min
Col. Temp. : 40 $^\circ\text{C}$
Detection : UV 254 nm
Injection Vol. : 0.5 μL

Analyte:

1. Ganciclovir
2. Aciclovir
3. Penciclovir

(100 mg/L each)