


Analysis of Antiviral drug nucleosides

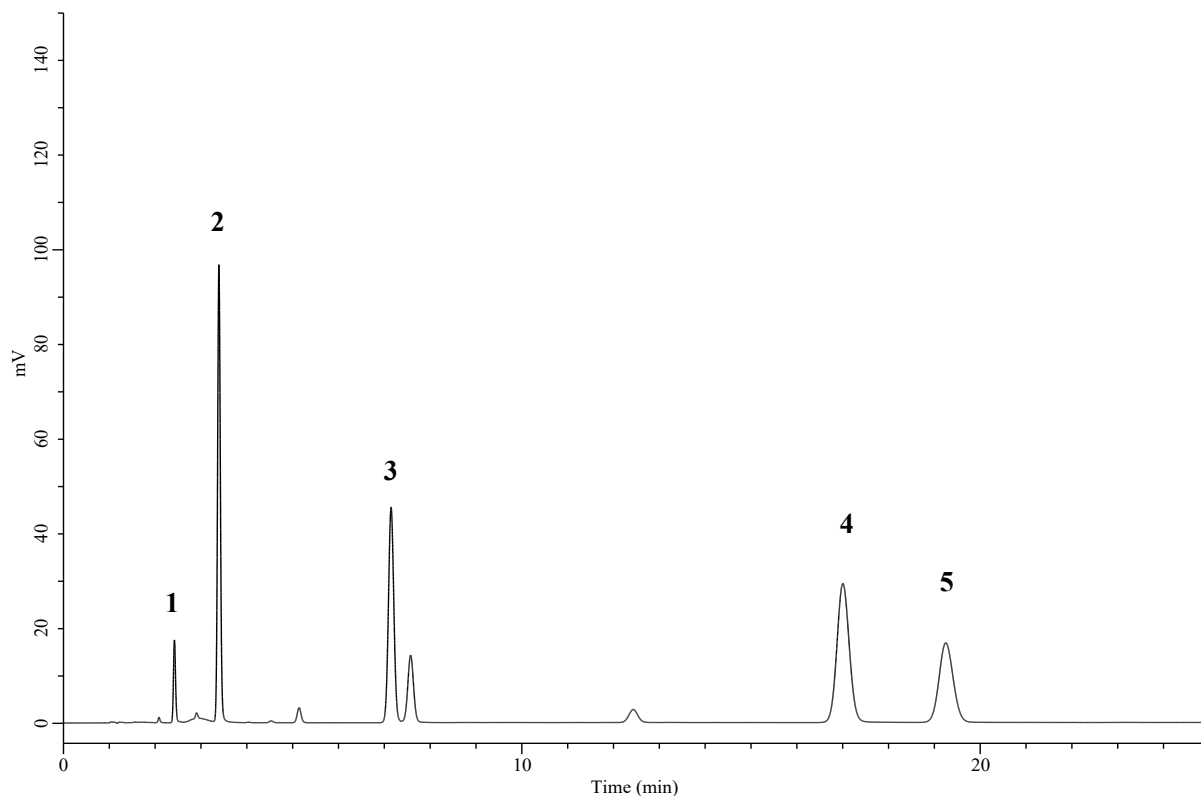

Conditions

Column : Inertsil Hybrid-C18 (GL Sciences Inc.)
(3.5 μm , 150 x 2.1 mm I.D.)
Column Cat. No. : 5020-76030
Eluent : 10 mM HCOONH_4 in H_2O
Flow Rate : 0.3 mL/min
Col. Temp. : 40 $^{\circ}\text{C}$
Detection : UV 240 nm
Injection Vol. : 1 μL

Analyte:

1. Ribavirin
2. Cytarabine
3. Gemcitabine
4. Lamivudine
5. Vidarabine

(100 mg/L each)