


Analysis of Antiviral drugs

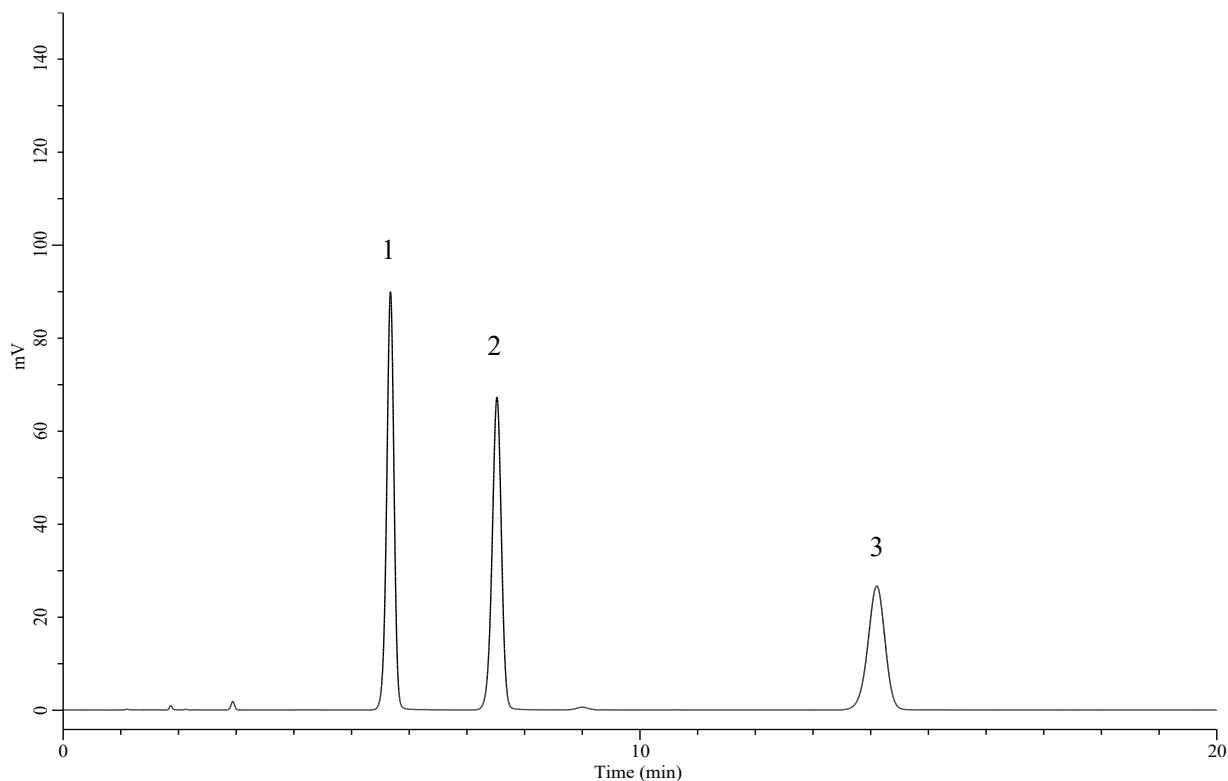

Conditions

System : Chromaster PLUS (HITACHI)
Column : Inertsil Hybrid-C18 (GL Sciences Inc.)
 (5 μ m, 150 x 2.1 mm I.D.)
Column Cat. No. : 5020-76004
Eluent : 10 mM HCOONH₄ in H₂O
Flow Rate : 0.3 mL/min
Col. Temp. : 40 °C
Detection : UV 254 nm
Injection Vol. : 1 μ L

Analyte:

1. Ganciclovir
2. Aciclovir
3. Penciclovir
(100 mg/L each)