


Analysis of Oligonucleotides

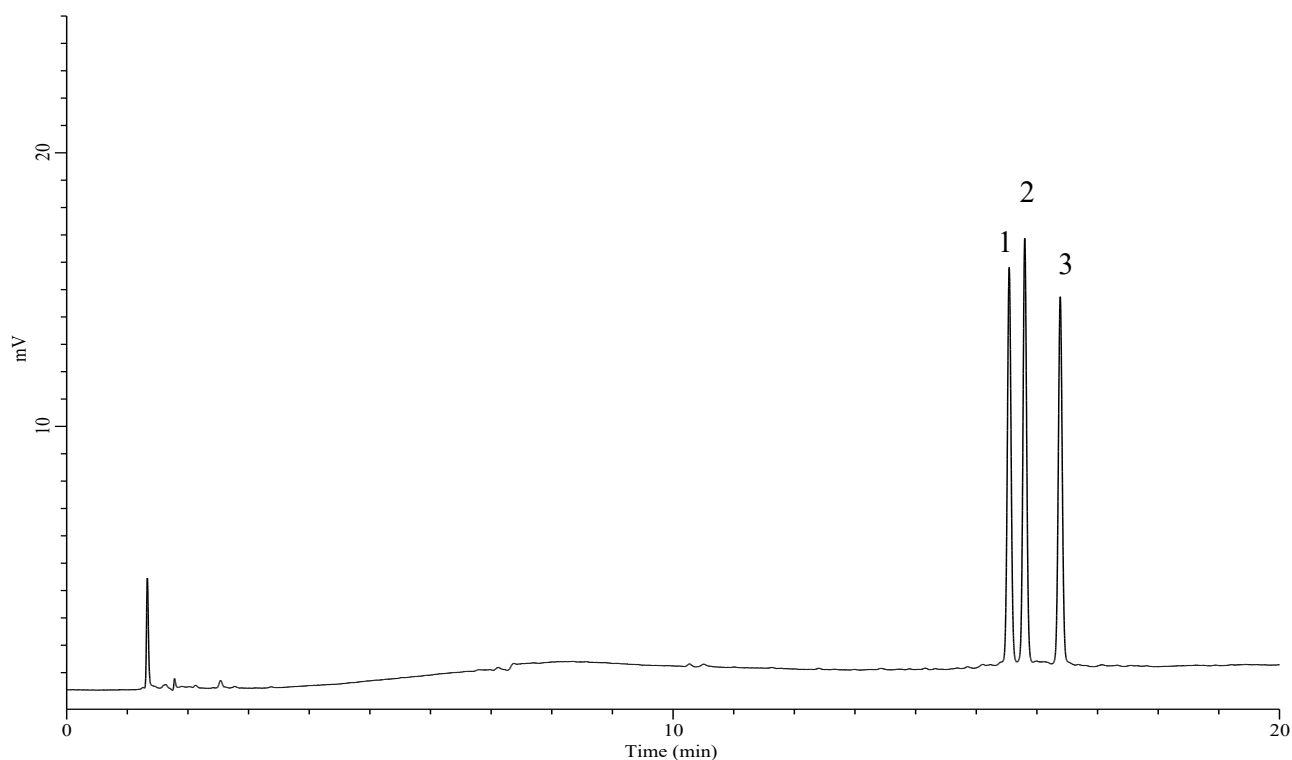

Conditions

Column : Inertsil Hybrid-C18 (GL Sciences Inc.)
(3.5 μm , 150 x 2.1 mm I.D.)

Column Cat. No. : 5020-76030

Eluent : A) 10 mM Triethylamine + 100 mM HFIP in ($\text{CH}_3\text{CN}/\text{H}_2\text{O}$ = 80/20, v/v)
B) 10 mM Triethylamine + 100 mM HFIP in H_2O

Time (min)	A (vol%)	B (vol%)
0	1	99
20	90	10

Flow Rate : 0.3 mL/min

Col. Temp. : 40 $^{\circ}\text{C}$

Detection : UV 254 nm

Injection Vol. : 1 μL

Analyte:

1. 5'-CATGACGTTTCCTGATGCT-3' 18 mer
2. 5'-CCATGACGTTTCCTGATGCT-3' 19 mer
3. 5'-TCCATGACGTTTCCTGATGCT-3' 20 mer
(10mg/L each)