


Analysis of Oligonucleotides

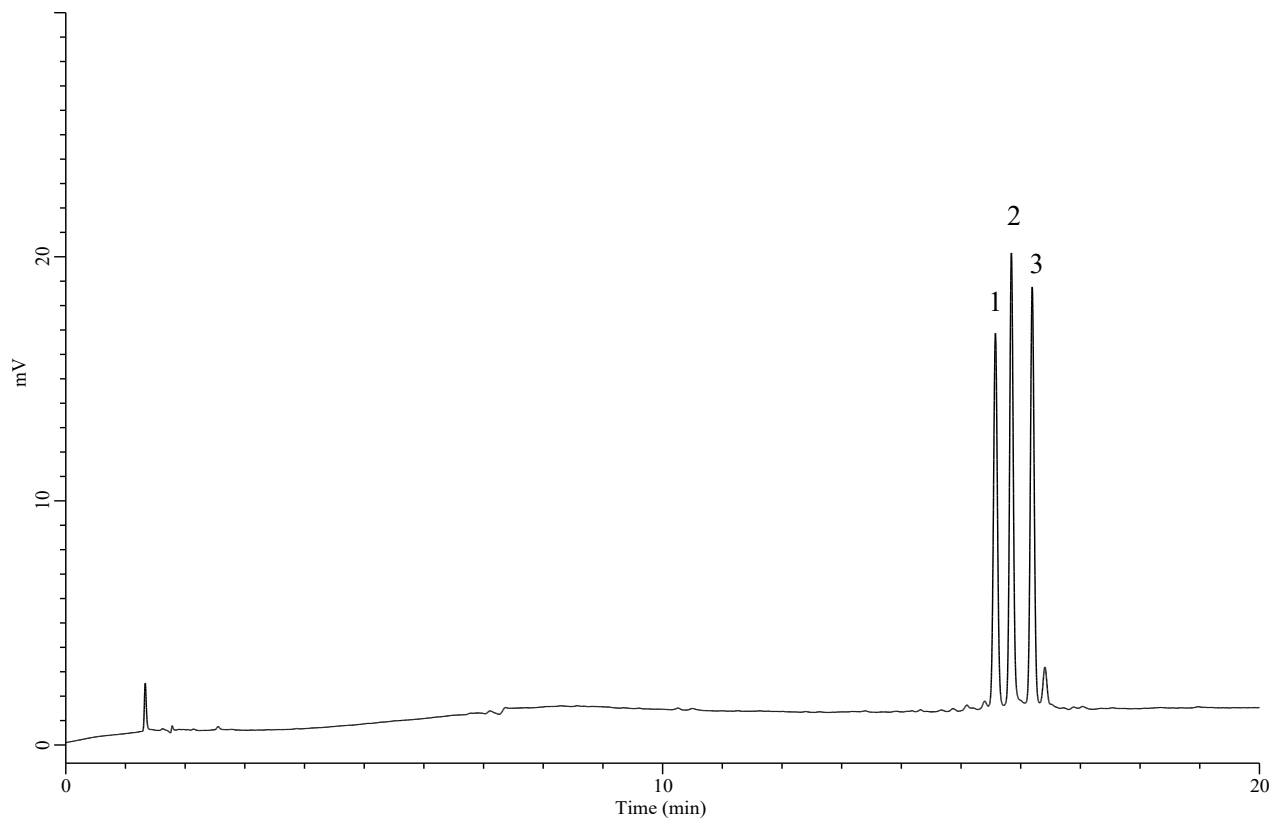

Conditions

Column : Inertsil Hybrid-C18 (GL Sciences Inc.)
(3.5 μm, 150 x 2.1 mm I.D.)
Column Cat. No. : 5020-76030
Eluent : A) 10 mM Triethylamine + 100 mM HFIP in (CH₃CN/H₂O = 80/20, v/v)
B) 10 mM Triethylamine + 100 mM HFIP in H₂O

Time (min)	A (vol%)	B (vol%)
0	1	99
20	90	10

Flow Rate : 0.3 mL/min
Col. Temp. : 40 °C
Detection : UV 254 nm
Injection Vol. : 1 μL

Analyte:

- 1. 5'- CATGACGTTTCCTGATGCT-3' 18 mer
 - 2. 5'- GCATGACGTTTCCTGATGCT-3' 19 mer
 - 3. 5'- TGCATGACGTTTCCTGATGCT-3' 20 mer
- (10 mg/L each)