


Analysis of Oligonucleotides

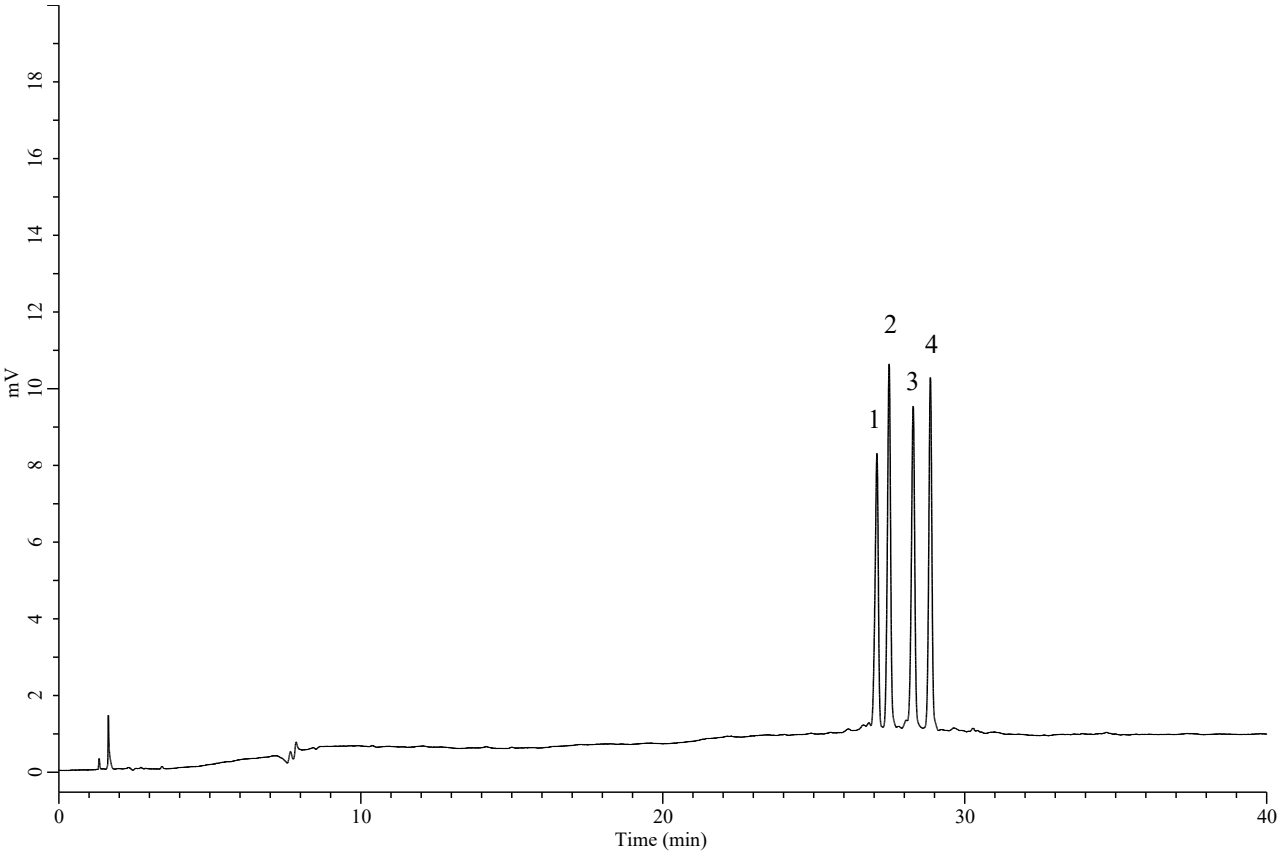

Conditions

Column : Inertsil Hybrid-C18 (GL Sciences Inc.)
(3.5 μ m, 150 x 2.1 mm I.D)
Column Cat. No. : 5020-76030
Eluent : A) 10 mM Triethylamine + 100 mM HFIP in (CH₃CN/H₂O =80/20, v/v)
B) 10 mM Triethylamine + 100 mM HFIP in H₂O

Time (min)	A (vol%)	B (vol%)
0	1	99
40	90	10

Flow Rate : 0.3 mL/min
Col. Temp. : 40 °C
Detection : UV 254 nm
Injection Vol. : 1 μ L

Analyte:

- 1. 5'- CATGACGTTTCCTGATGCT-3' 18 mer
- 2. 5'- TATGACGTTTCCTGATGCT-3' 18 mer
- 3. 5'- GCTAGGACGTTTCCTGATGCT-3' 20 mer
- 4. 5'- GCTATGACGTTTCCTGATGCT-3' 20 mer

(10 mg/L each)