


Analysis of Steroids

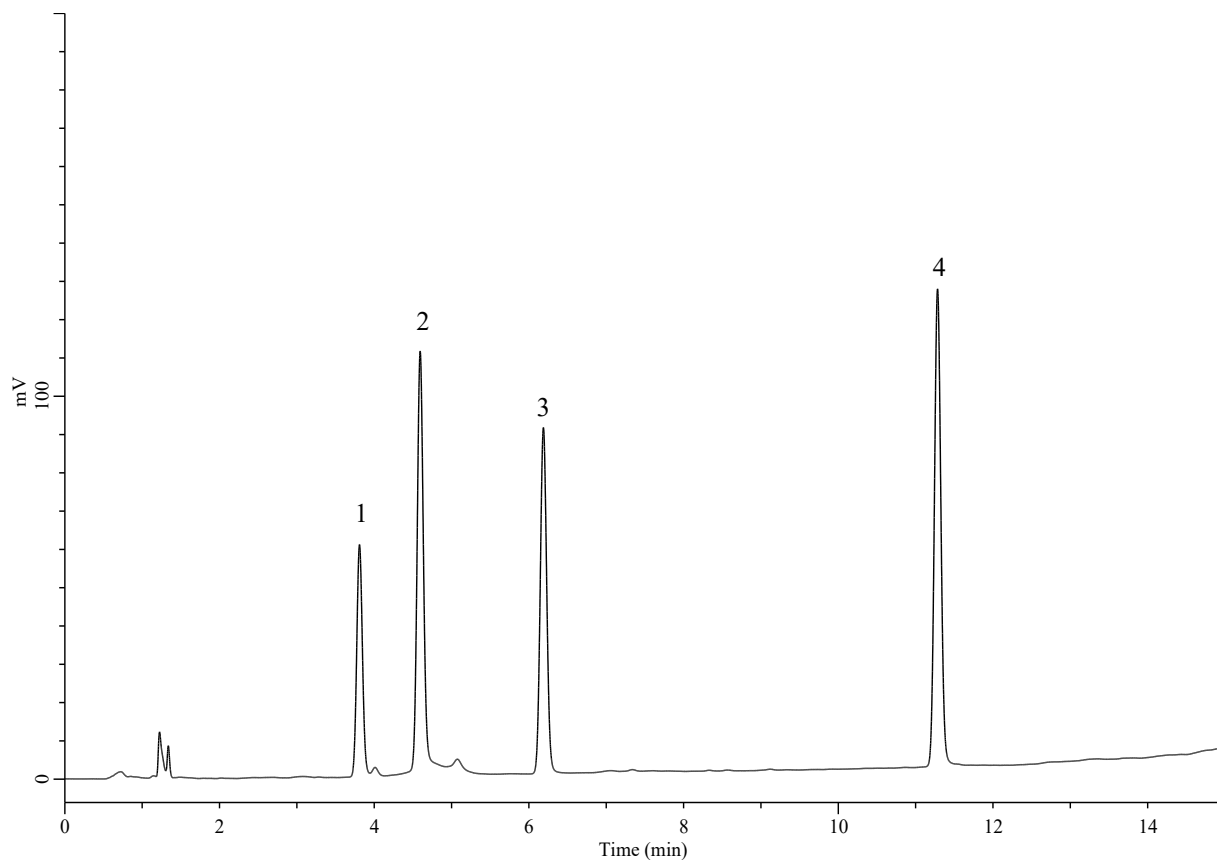

Conditions

Column : Inertsil Hybrid-C18 (GL Sciences Inc.)
(3.5 μm , 150 x 2.1 mm I.D.)

Column Cat. No. : 5020-70594

Eluent : A) CH_3OH
B) H_2O

Time (min)	A (vol%)	B (vol%)
0	50	50
40	90	10

Flow Rate : 0.3 mL/min

Col. Temp. : 40 $^{\circ}\text{C}$

Detection : UV 240 nm

Injection Vol. : 1 μL

Analyte:

1. Prednisolone
2. Prednisolone
3. Corticosterone
4. Progesterone

(100 mg/L each)