


Analysis of Heparin sodium

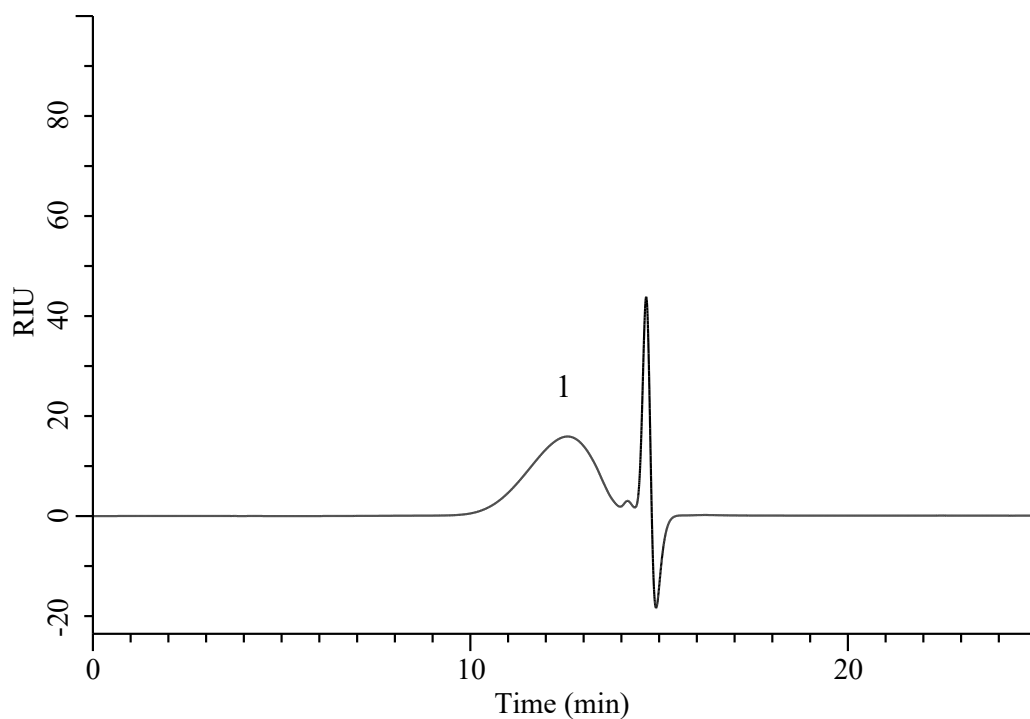

Conditions

System : Chromaster PLUS (HITACHI)
Column : Inertsil WP300 Diol (GL Sciences Inc.)
 (5 μ m, 250 x 7.6 mm I.D.)
Column Cat. No. : 5020-05988
Eluent : 0.2 M $\text{CH}_3\text{COONH}_4$ in H_2O
Flow Rate : 0.6 mL/min
Col. Temp. : 40 $^\circ\text{C}$
Detection : RI (Chromaster 5450 RI)
Injection Vol. : 20 μL
Sample : Standard in H_2O

Analyte:

1. Heparin sodium,
from porcine intestinal mucosa 10 g/L
(mol wt 4000-6000 Da)