


Analysis of Dotinurad

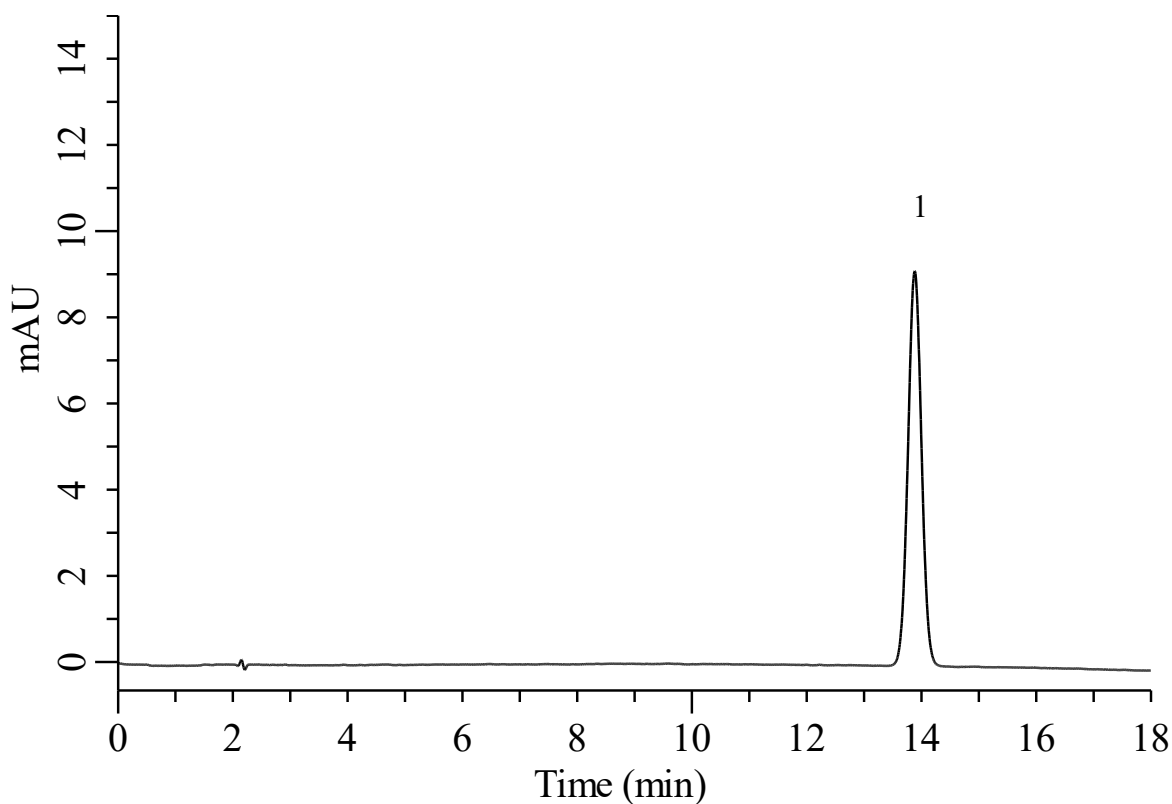

Conditions

System : Primaide (HITACHI)
Column : InertSustain C18 (GL Sciences Inc.)
 (5 μ m, 250 x 4.6 mm I.D.)
Column Cat. No. : 5020-07346
Eluent : A: CH₃CN
 B: 0.1% HCOOH in H₂O
 A/B = 40/60, v/v
Flow Rate : 1.0 mL/min
Col. Temp. : 40 °C
Detection : UV 250 nm (Primaide 1435 DAD)
Injection Vol. : 10 μ L
Sample : Standard in Eluent

Analyte:

1. Dotinurad 10 mg/L