


Analysis of Dotinurad

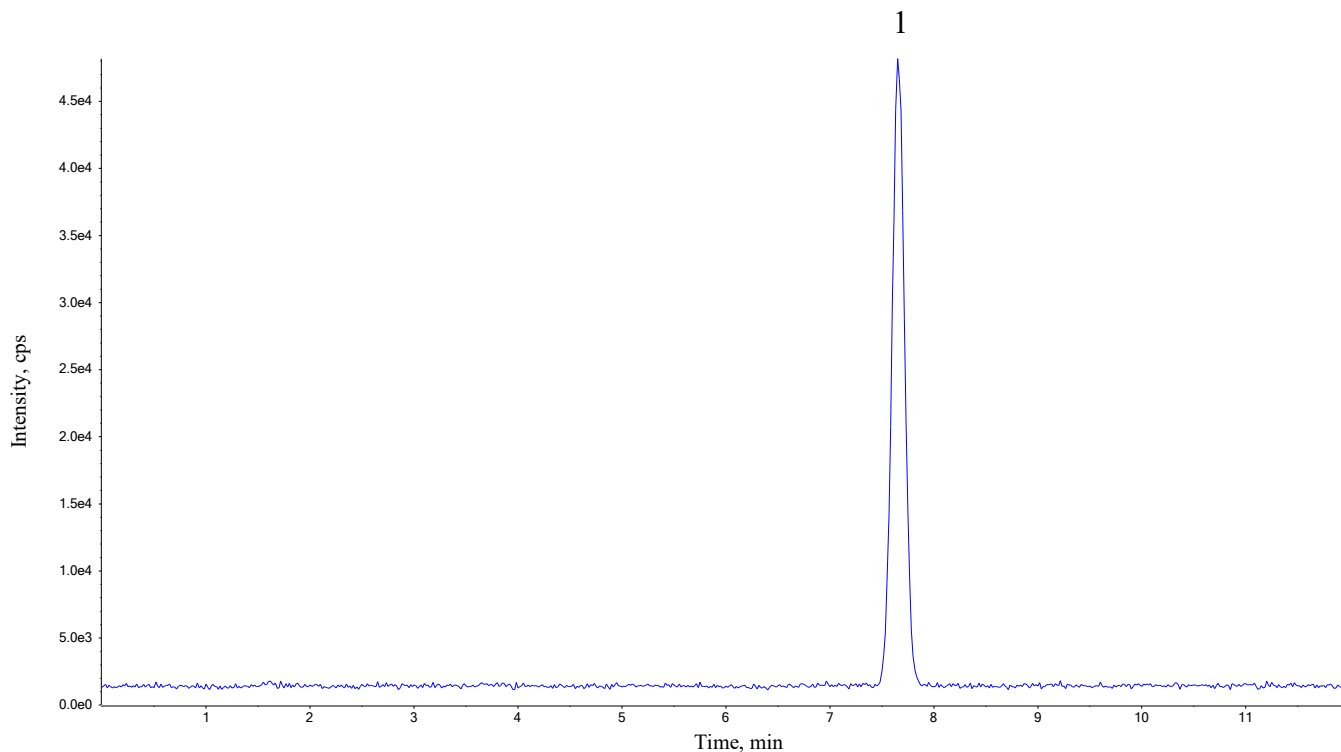

Conditions

System	: ExionLC & Triple Quad 4500 (Sciex)
Column	: Inertsil Hybrid-C18 (GL Sciences Inc.) (3.5 μm, 150 x 2.1 mm I.D.)
Column Cat. No.	: 5020-76030
Eluent	: A: CH ₃ CN B: 0.1% HCOOH in H ₂ O A/B = 40/60, v/v
Flow Rate	: 0.2 mL/min
Col. Temp.	: 40 °C
Injection Vol.	: 2 μL
Sample	: Standard in Eluent
Detection	: MS/MS (ESI, Negative, MRM) CUR CAD IS TEM GS1 GS2 20 7 -4500 400 50 50

Analyte:

1.	Dotinurad	10 µg/L	Q1 / Q3 357.6 / 161.7
----	-----------	---------	--------------------------