


Analysis of Antiviral Drugs

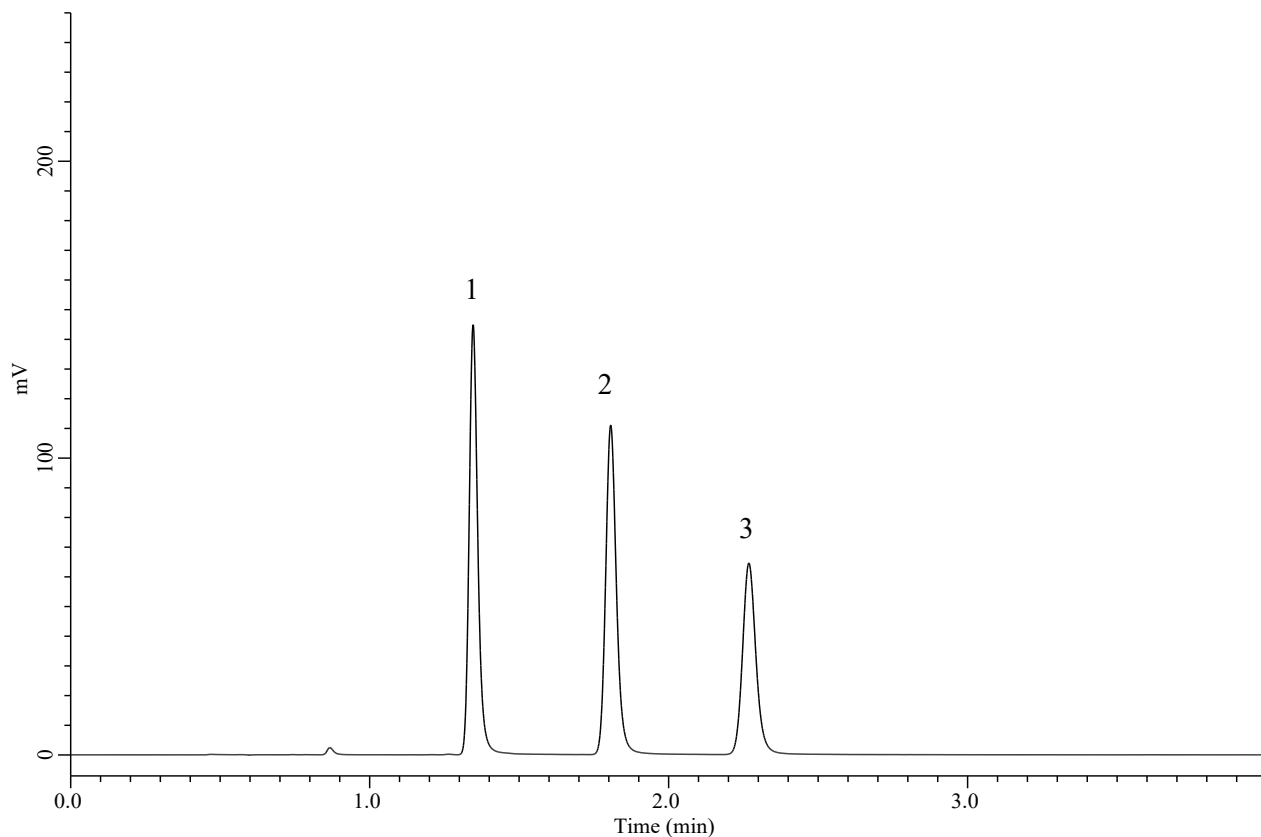

Conditions

Column : InertCore PFP
 (2.6 μ m, 100 x 2.1 mm I.D)
Column Cat. No. : 5020-17559
Eluent : 0.1% HCOOH + 10 mM HCOONH₄ in H₂O
Flow Rate : 0.4 mL/min
Col. Temp. : 40 °C
Detection : UV 254 nm
Injection Vol. : 0.5 μ L

Analyte

1. Ganciclovir
 2. Aciclovir
 3. Penciclovir
 (100 mg/L of each)