


Analysis of Thymol and Carvacrol

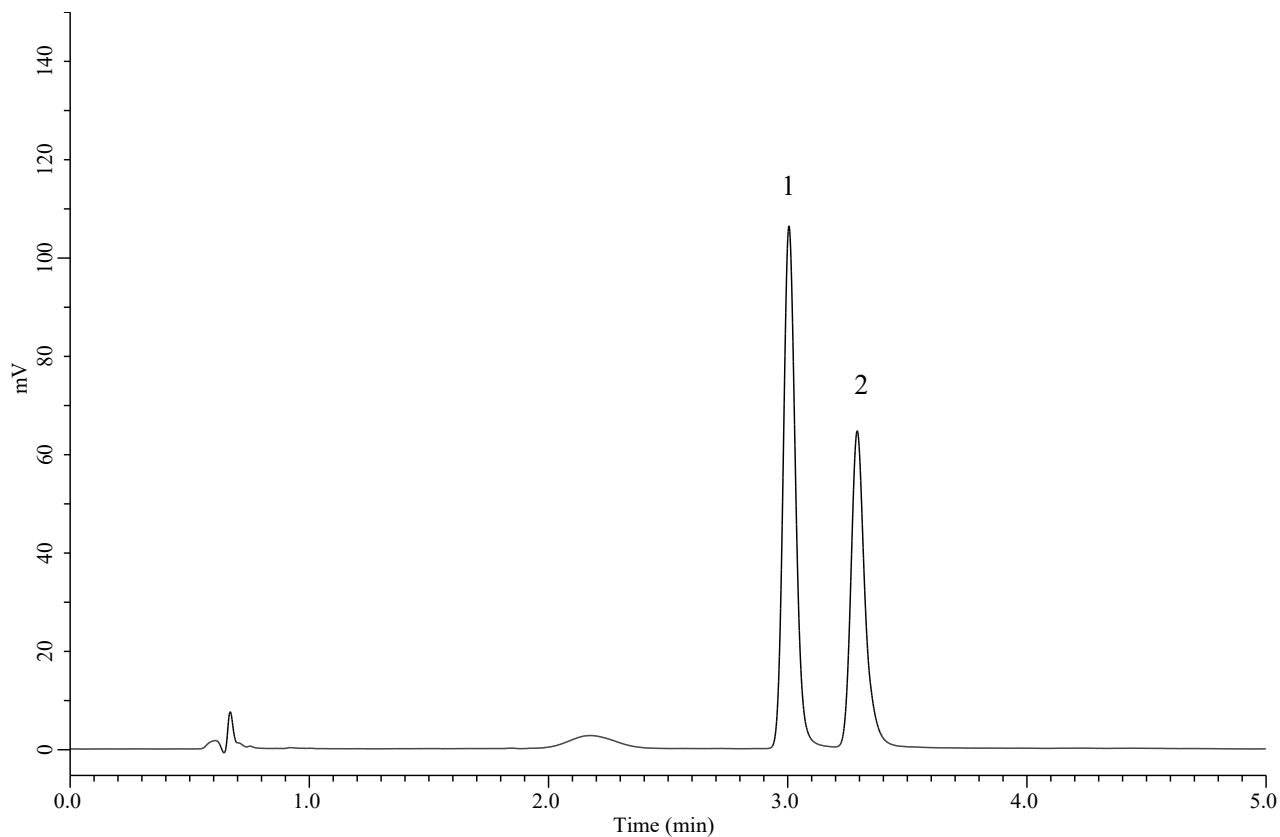

Conditions

Column : InertCore PFP
 (2.6 μ m, 100 $\times$ 2.1 mm I.D)
Column Cat. No. : 5020-17559
Eluent : A) CH₃CN
 B) 0.1% HCOOH in H₂O
 A/B = 35/65, v/v
Flow Rate : 0.3 mL/min
Col. Temp. : 40 °C
Detection : UV 275 nm
Injection Vol. : 1 μ L

Analyte

1. Thymol
 2. Carvacrol
 (100 mg/L of each)