


Analysis of Vitamin E

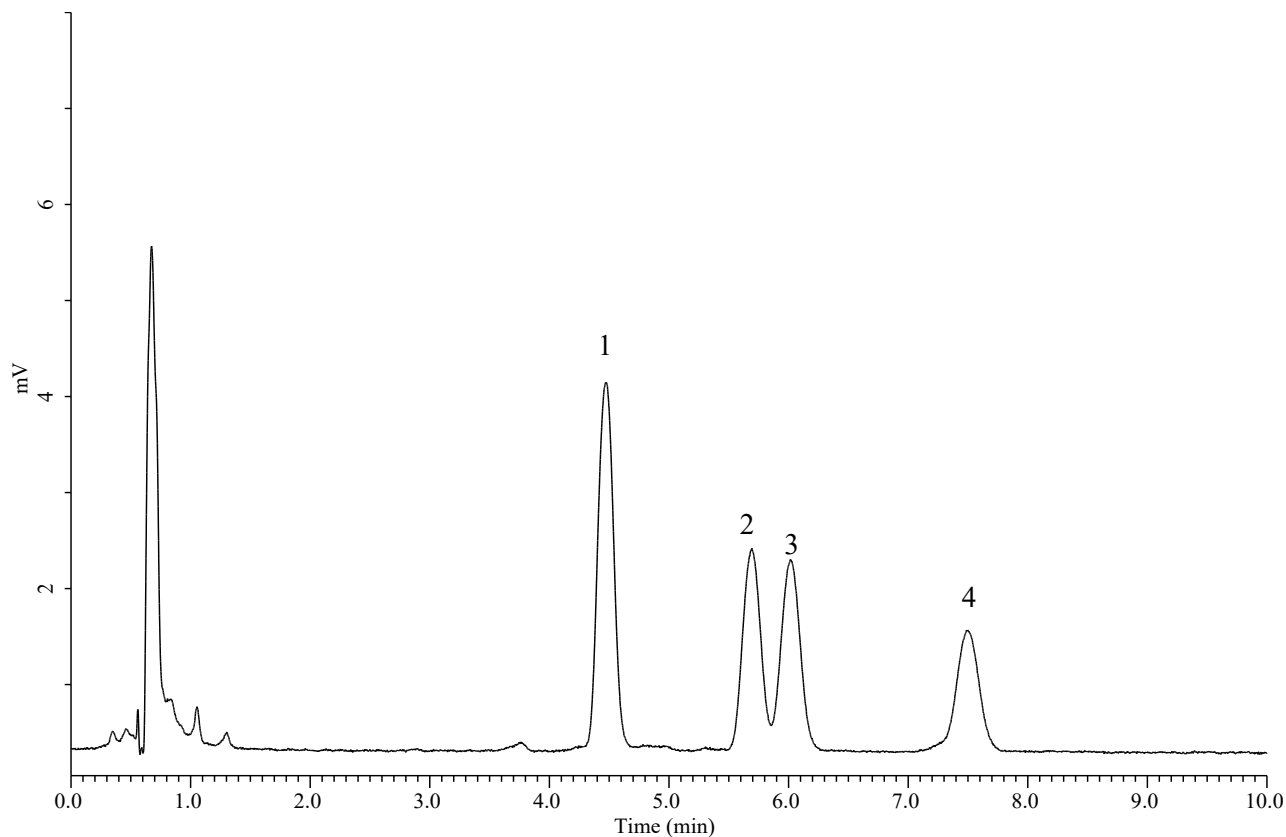

Conditions

Column : InertCore PFP
 (2.6μm, 150 × 2.1 mm I.D)
Column Cat. No. : 5020-17560
Eluent : A) CH₃OH
 B) H₂O
 A/B = 85/15, v/v
Flow Rate : 0.3 mL/min
Col. Temp. : 40 °C
Detection : UV 290 nm
Injection Vol. : 2 μL

Analyte

1. α-Tocopherol
 2. β-Tocopherol
 3. γ-Tocopherol
 4. δ-Tocopherol
 (100 mg/L of each)