


Analysis of *cis*, *trans*-Citral

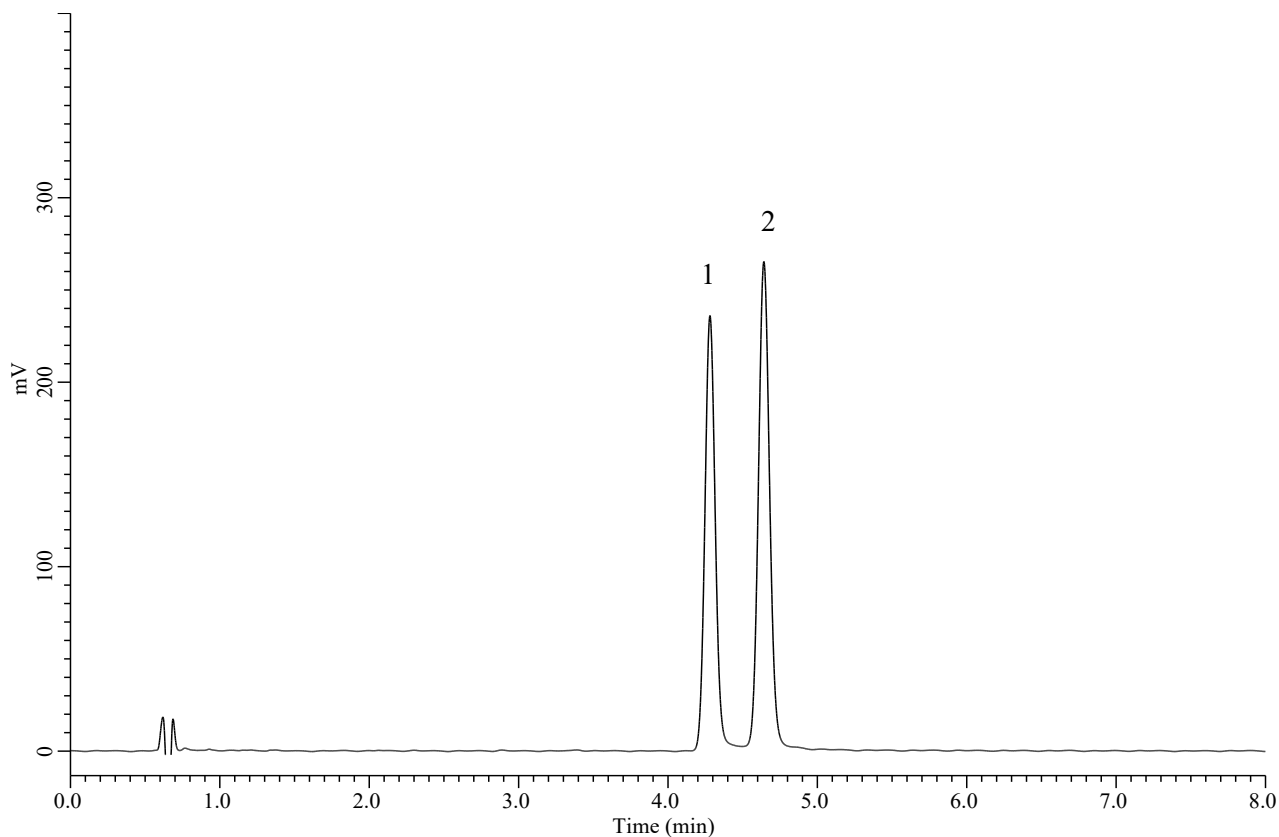

Conditions

Column : InertCore PFP
 (2.6 μ m, 100 $\times$ 2.1 mm I.D)
Column Cat. No. : 5020-17559
Eluent : A) CH₃CN
 B) 0.1% HCOOH in H₂O
 A/B = 35/65, v/v
Flow Rate : 0.3 mL/min
Col. Temp. : 40 °C
Detection : UV 240 nm
Injection Vol. : 1 μ L

Analyte

1. *cis*-Citral (Neral)
 2. *trans*-Citral (Geranial)
 (500 mg/L of each)