


Analysis of Antiviral Drugs

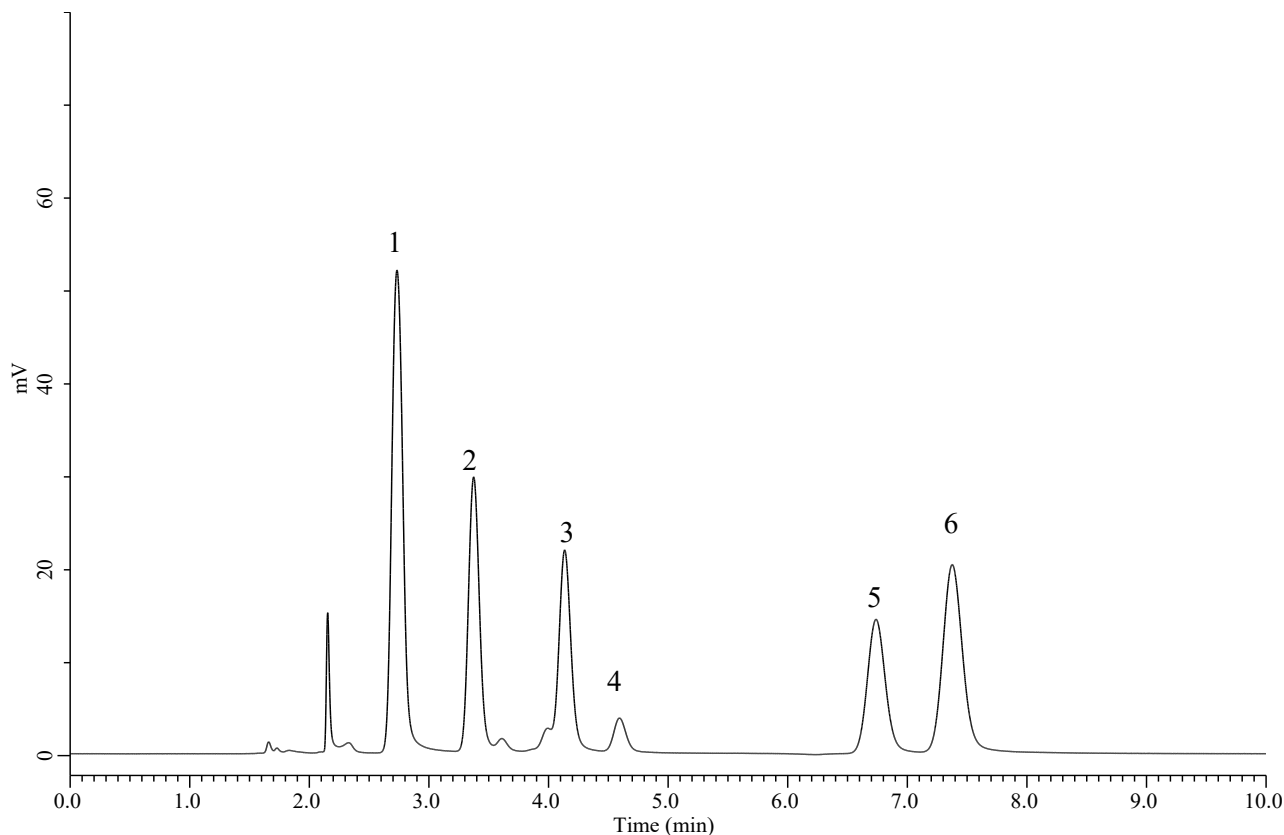

Conditions

Column : Inertsil Hybrid-Z HILIC
 (3.5 μ m, 150 $\times$ 2.1 mm I.D)
Column Cat. No. : 5020-76196
Eluent : A) CH₃CN
 B) 10 mM HCOONH₄ in H₂O
 A/B = 90/10, v/v
Flow Rate : 0.3 mL/min
Col. Temp. : 40 °C
Detection : UV 240 nm
Injection Vol. : 1 μ L

Analyte

1. Lamivudine
 2. Gemcitabine
 3. Vidarabine
 4. Ribavirin
 5. Cytarabine
 (100 mg/L of each)