


Analysis of Sugar Alcohols

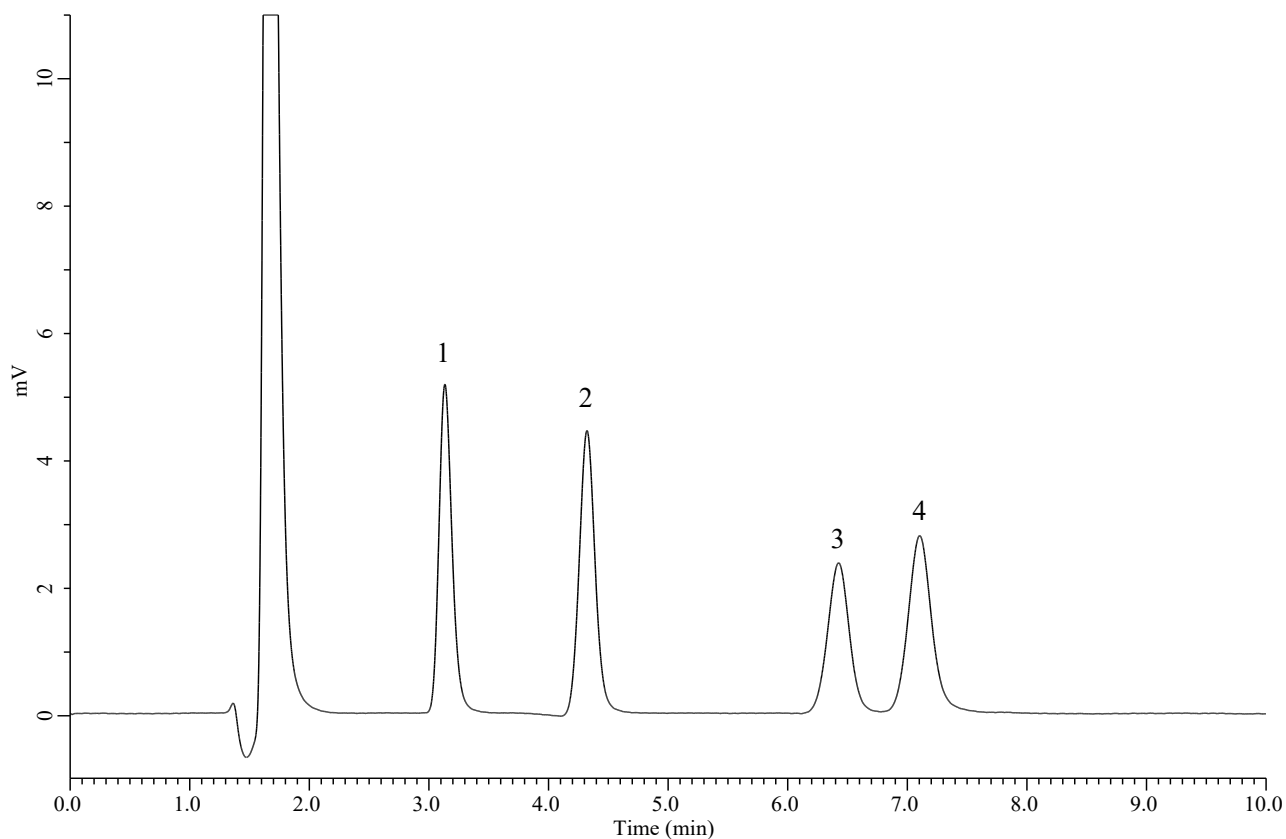

Conditions

Column : Inertsil Hybrid-Z HILIC
(5 μm , 150 $\times$ 2.1 mm I.D)
Column Cat. No. : 5020-76173
Eluent : A) CH_3CN
B) 10 mM HCOONH_4 in H_2O
A/B = 90/10, v/v
Flow Rate : 0.3 mL/min
Col. Temp. : 40 $^\circ\text{C}$
Detection : RI
Injection Vol. : 4 μL

Analyte

1. Erythritol
2. Xylitol
3. Sorbitol
4. Mannitol
(5000 mg/L of each)