


Analysis of Organic solvents (Benzene, Toluene, Ethylbenzene, Xylene)

Specification

One GC Channel

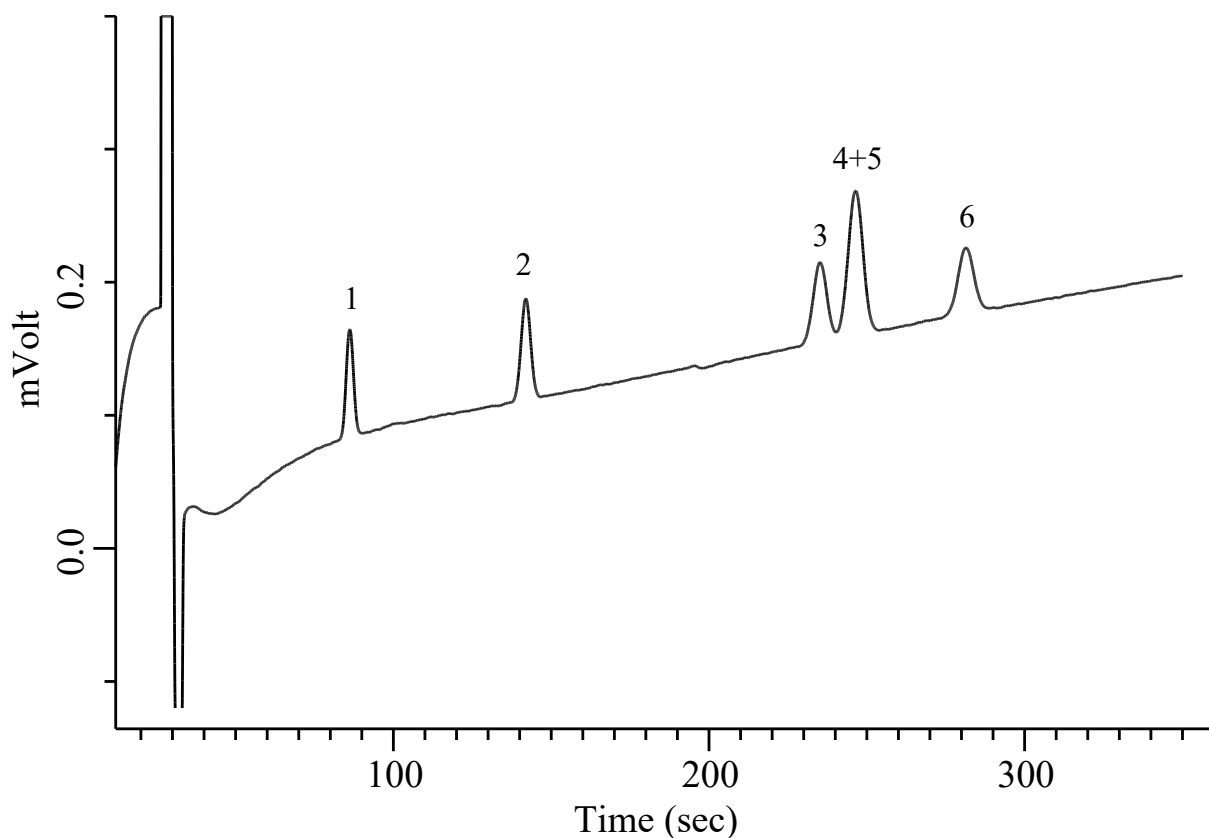

Conditions

System : Micro GC
Column : CP-Volamine 15 m
Column Temp. : 100 °C
Carrier Gas : He 190 kPa
Inject Time : 400 msec

Analyte:

1. Benzene	10 ppm (v/v)
2. Toluene	10 ppm (v/v)
3. Ethylbenzene	10 ppm (v/v)
4. <i>m</i> -Xylene	10 ppm (v/v)
5. <i>p</i> -Xylene	10 ppm (v/v)
6. <i>o</i> -Xylene	10 ppm (v/v)

Balance Gas N₂